


**Nekteck**

# Instruction Manual

## Air Compression Full Zipper Leg Massager

Model No.: FE-7204B

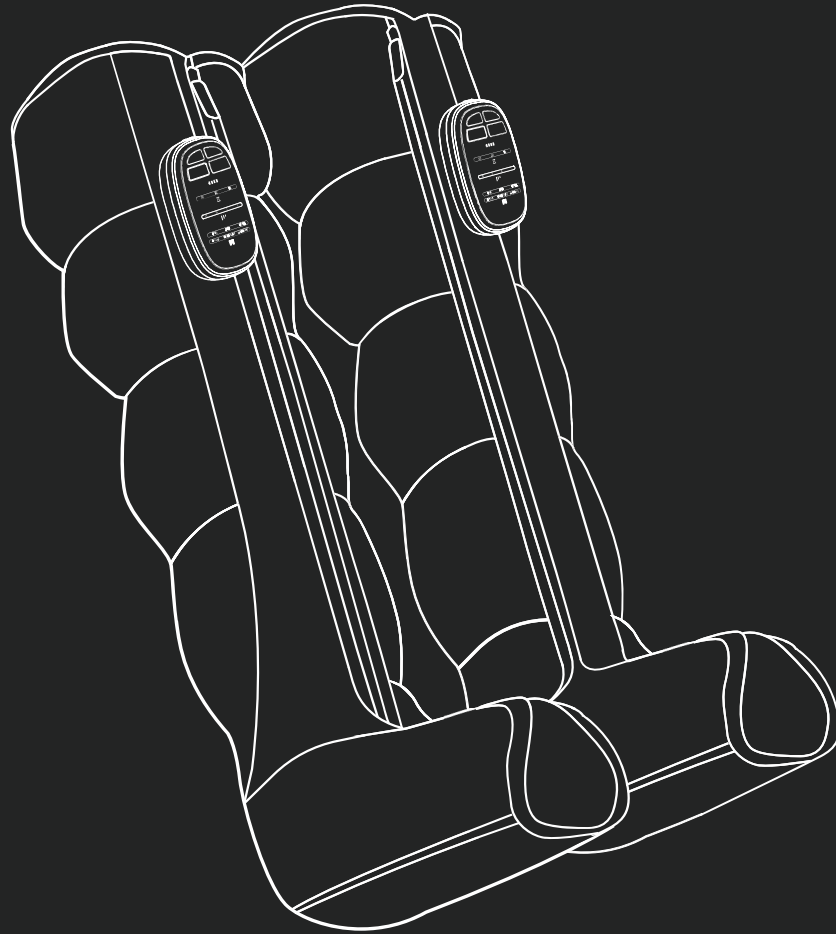

**Questions or Concerns?**

Contact us at [support@nekteck.com](mailto:support@nekteck.com)



**TABLE OF CONTENTS**

Important Safeguards .....	01
Getting To Know Your Air Compression Full Zipper Leg Massager .....	03
Before First Use .....	05
Operating Instructions .....	06
FAQS .....	08
Care & Maintenance .....	10
Troubleshooting .....	11
FCC Statement .....	12
Instruction For Reporting Adverse Event .....	13
Submitting Adverse Event Reports To FDA .....	13
Warranty Information .....	14
Customer Service .....	16

**TABLE DES MATIÈRES**

Mesures De Sécurité Importantes .....	18
Découvrez votre Appareil de Pressothérapie Pour Jambes à Fermeture Éclair .....	20
Avant La Première Utilisation .....	22
Fonctionnement de l'appareil.....	23
FAQ .....	25
Entretien Et Maintenance .....	27
Dépannage .....	28
Déclaration de conformité FCC .....	28
Instructions Pour Signaler Un Événement Indésirable .....	30
Soumission Des Rapports D'événements Indésirables À La FDA .....	30
Informations Sur La Garantie .....	31
Service Clientèle .....	33

# IMPORTANT SAFEGUARDS

## PRECAUTIONS TO ENSURE SAFE AND RELIABLE OPERATION:

- Do not insert any foreign objects into the openings.
- Do not attempt to disassemble the unit.
- Do not crumple or crush the massager by any means.
- Do not place the unit in direct sunlight nor store in dusty places or with high humidity.
- Ensure that children do not play the unit unsupervised; some parts may be swallowed.
- Do not use the device near strong electrical or electromagnetic fields generated by cell phones or other devices, as they may cause or become a source of interference.
- Any modification to the unit shall not be allowed.
- The product (Air Compression Full Zipper Leg Massager) is for adult use only.

## CONTRAINDICATIONS:

This Air Compression Full Zipper Leg Massager is not intended for use if you have any of the following medical conditions. We advise you to consult with a physician first to confirm if the product is clinically safe for you to use.

- Presumptive evidence of congestive heart failure as increased blood return to the heart could be detrimental.
- Pre-existing deep vein thrombosis (DVT) or blood clots.
- Existing hemorrhagic disease or coagulation problems.
- Thrombophlebitis
- Pulmonary embolism

This product should not be utilized by individuals with edema, ischemic vascular disease, or severe peripheral neuropathy.

## SAFETY PRECAUTION:

- Do not attempt to repair or dismantle the product by yourself or a third party. This will void the warranty.
- Do not use the product in a damp place or where it may get wet. If the unit is wet, do not use it.
- Stop using the product immediately if you feel uncomfortable or any abnormalities occur.

- Do not use this product in a high-temperature environment.
- Do not place the unit in a place where it may be in contact with flammable items, explosive items, or chemicals.
- Do not touch the unit with any sharp object or risk damage.
- Do not use this product if any of the below apply to you: Skin is infected, swollen, burnt or otherwise damaged. The product should be used on healthy skin.

Please do not use the product when you have the following medical conditions:

- You have a pacemaker or any other medical implant
- You are suffering from heart disease.
- You are suffering from diabetes.
- Your bones are fragile or injured.
- You are suffering from cancer.
- You are receiving treatment or your blood pressure is not normal.
- You are suffering from osteoporosis.
- You are currently pregnant or having menstruation.
- Children or adolescents who are developing physically.




**SAVE THESE INSTRUCTIONS  
HOUSEHOLD USE ONLY**


# GETTING TO KNOW YOUR AIR COMPRESSION FULL ZIPPER LEG MASSAGER

Air Compression Full Zipper Leg Massager, is intended for the temporary relief of minor muscle aches and pains and for the temporary increase in circulation to the treated foot/calf and thigh in people who are in good health. The Air Compression Full Zipper Leg Massager simulates kneading and stroking of tissues by using air-inflatable boots.

## 1) PRODUCT DIAGRAM

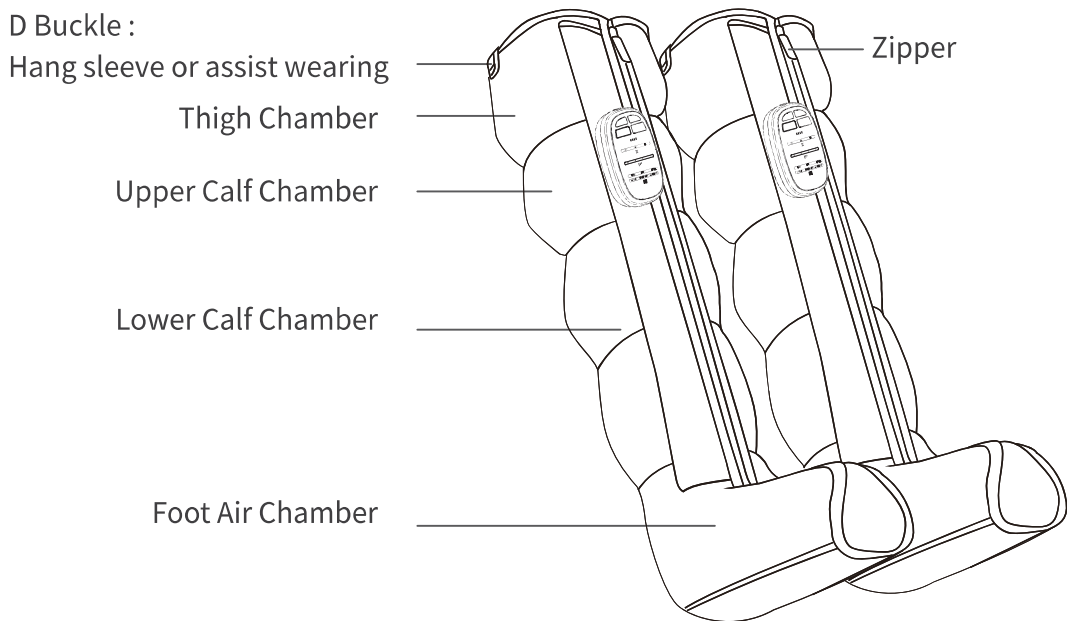


## 2) CONTROL PANEL

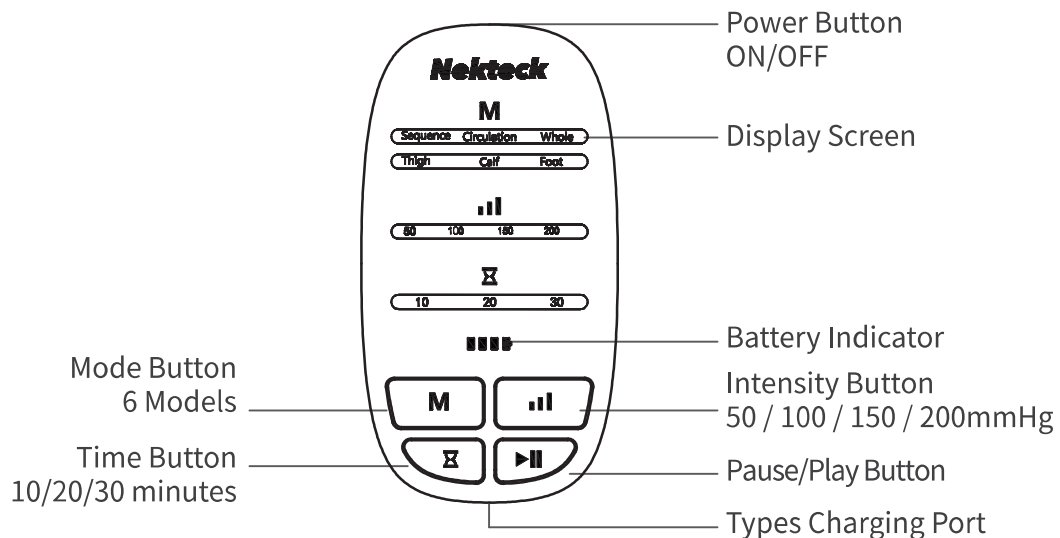

**3) PACKAGE CONTENTS**

- 2 x Air Compression Full Zipper Leg Massager
- 2 x USB Charging Cable
- 1 x Portable Bag
- 1 x User Manual

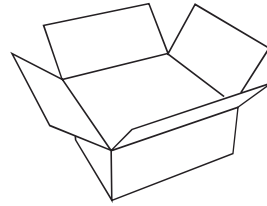**4) SPECIFICATIONS**

<b>Model No.</b>	FE-7204B
<b>DC output</b>	5V 2A
<b>Rated Power</b>	10W
<b>battery capacity</b>	3000mAh
<b>Unit size (LxWxH)</b>	15.75 x 7.87 x 8.66 inch
<b>Weight</b>	7.15 Lbs
<b>Leg Boots Size</b>	Inseam: 33in; Max Thigh: 30in; Max Foot Width: 18in
<b>Website</b>	<a href="http://www.nekteck.com">www.nekteck.com</a>
<b>Email</b>	<a href="mailto:support@nekteck.com">support@nekteck.com</a>

# BEFORE FIRST USE

## USING THE MASSAGER

1. Remove the Air Compression Full Zipper Leg Massager and its accessories from the packaging box. Ensure all items and accessories are complete.
2. Check that the fabric is in good condition.
3. Sit in your most comfortable posture, place your feet and legs into the massager and press the control to start the Air Compression massage. Adjust the levels to your preferred comfort.
4. After turning it on, be sure to check whether the device is working properly, and then select your preferred massage mode and intensity. Adjust the levels to your preferred comfort.
5. There may be signs of use on the skin of the boots. If you are concerned about signs of use, wear socks or pants before using the massager. Please be aware that this device can leave temporary marks on your feet, calves, knees, and thighs.






**Note:** Do not use on zippered clothes, pants or metal objects. Otherwise the product may be damaged.

## TIPS:

Ensure the zipper is completely zipped up before wearing. Never attempt to use the system with the zipper partially or fully unzipped.

- \* You can use the D-buckle on top of the sleeve to help put it on quickly.

# I OPERATING INSTRUCTIONS

BUTTON	FUNCTION
 POWER BUTTON	<p>Press this button for 2 seconds to turn on the device. The default setting will be: Mode 1(Sequence Mode), 50 intensity(the lowest), 10 minutes timing.</p> <p>Press it again for 2 seconds to turn off the device.</p>
 MODE BUTTON	<p>Press mode button to change the massage mode. 6 modes optional,see the difference in FAQs A2 session.</p> <p>Mode 1 : Sequence Mode 2 : Circulation Mode 3 : Whole Mode 4 : Thigh Mode 5 : Calf Mode 6 : Foot</p>
 INTENSITY BUTTON	<p>Press intensity button to adjust air pressure intensity, 4 levels available (50 / 100 / 150 / 200mmHg).</p> <p>* We suggest starting with low intensity for your first use and gradually increasing it as you get used to it.</p>
 TIMER BUTTON	<p>Press the time button to change between 3 treatment session times, 10 mins, 20 mins or 30 mins. The time can be adjusted anytime while the session is running.</p>
 PAUSE BUTTON	<p>Press the Pause button once to pause the product. In the pause state, all function buttons except Pause and Power will be inactive. Press Pause again to resume operation.</p> <p>* The pause function allows you to interrupt the massage program at any time to adjust intensity or handle other matters. This provides a more flexible and personalized experience, giving users greater control over the massage's effect and rhythm.</p>

**BATTERY INDICATOR:**

1. During use, the battery icon will progressively turn off as the battery power decreases.
2. When charging, the battery icon will light up sequentially. Charging is complete when the battery is full, and the indicator light will show a steady white light. Once fully charged, unplug the charging cable, and the battery icon light will turn off.

**MEMORY FUNCTION:**

The device will remember the current working mode status (mode, intensity, timer and gear) after starting the machine. It will become the default setting when restarting the machine the next time.

**PRESSURE SENSOR:**

Precise control of pressure intensity for each level adapts to the different demands of various users, enhancing their experience.






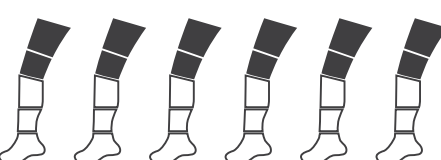
# FAQS

## Q1 How does this product work?

**A1** This is an air compression device with 5 air chambers (foot, lower calf, upper calf, lower thigh, and upper thigh) that inflate and deflate in a sophisticated sequence to simulate hand kneading and stroking. It maximizes blood flow throughout the body, helping to reduce muscle soreness and pain in the treated areas, and improve circulation.

## Q2 How many massage modes, and what's the difference?

**A2** There are 6 massage modes. Please check the differences below.

Sequence	
Circulation	
Whole	
Foot	
Calf	
Thigh	

**Q3 What should I do if I feel the massage strength is too weak or too strong?**

**A3** There are 4 levels of massage intensity available. Please select the intensity that suits you. If you still feel it's too strong at the lowest setting, please stop using the device and consult your doctor before your next use.

**Q4 Under what circumstances should I not use it?**

**A4** During pregnancy or immediately after giving birth. Patients with severe blood circulation disorders are not advisable to use the product under the above circumstances. In addition, if you have a fever or any abnormality in your body, be sure to consult a doctor before use.

**Q5 How many times can I use it in a day?**

**A5** The massager will automatically shut down after 30 minutes of continuous use. We do not recommend using a massager for more than 40 minutes in a day. Prolonged use will over-stimulate muscles and nerves and may have adverse effects.

**Q6 Why does the control unit make sound?**

**A6** The sound comes from the working air pump in the control unit, which continuously provides air to the airbags in the sleeves. It's a low sound, and that's a normal phenomenon.

**Q7 When should I use this product?**

**A7** This product promotes blood flow, effectively removing muscle soreness and fatigue. You can use it anytime you want: before a workout to warm up your muscles, after exercise for better and faster recovery, or when you get home after a long day at work.

**Note:**

- Expected service life: Approximately two years at 2 uses per day.  
Do not service or perform maintenance on the device while it is in use.  
Maintenance should only be performed by service personnel.
- Service and maintenance of your device may require repair, parts, and technical support. As noted in the beginning of the user manual, users should contact the manufacturer or the representative of the manufacturer for using or maintaining the product.

# CARE & MAINTENANCE

## USE ENVIRONMENT

	Use environment	Storage and transportation
Temperature	5°C~40°C	-10°C+60°C
Humidity range	5%~93%RH	5%~90%RH
Barometric range	86kPa~105kPa	86kPa~105kPa

## SECURITY AND MAINTENANCE

- Never use thinner, benzene, or alcohol to clean (as these may cause malfunction, damage to parts, or discoloration).
- Avoid dropping, slamming, or throwing the massager.
- Do not expose the unit directly under the sunlight.

## WET THE CLOTH AND WRING IT OUT

If you are concerned about dust, use a neutral kitchen detergent diluted to 5% with warm or clean water. Dampen a soft cloth with this solution and then wipe the unit clean.

## CLEAN

Make sure the power is turned off. Gently wipe the inside and outside of the leg massager, avoiding excessive force.

## NOTE:

If rubbish or hair gets caught in the Velcro, remove it with a toothpick or similar tool.

## STORAGE

Make sure the product is turned off and the gas is exhausted, carefully fold the air pressure leg massager and put it into the storage bag. Make sure the massager is not bent or deformed, which may affect next use. Close the storage bag for easy carrying and storage.

## SAFEKEEPING

Put the unit in a cool, dry place and avoid direct sunlight or high-temperature places.

- Keep out of the reach of children.
- Do not store other objects on the surface of the controller or massager to avoid damaging the device's surface and to avoid discoloration of the product.

# I TROUBLESHOOTING

Before requesting repairs, please confirm the following contents, and if errors persist after inspection, please contact [support@nekteck.com](mailto:support@nekteck.com).

CONDITIONS	TROUBLESHOOTING
The unit does not turn on	Check if the device is fully charged
Stop working halfway	<ol style="list-style-type: none"> <li>1. Check whether the device has insufficient power.</li> <li>2. The device has a scheduled shutdown function.</li> <li>3. Check if the pause button is pressed. (When paused, all functions are temporarily disabled, and button commands are invalid. Briefly press the pause button again to resume.)</li> </ol>
Strength intensity is too weak	<ol style="list-style-type: none"> <li>1. Adjust the massager intensity setting to increase the message intensity.</li> <li>2. Check whether the device has insufficient power and charge it.</li> </ol>
Strength intensity is too strong	Adjust the massager intensity setting to reduce the message intensity.

## FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

**NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

**FCC Caution:** Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

### FCC RADIATION EXPOSURE STATEMENT

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. End users must follow the specific operating instructions for satisfying RF exposure compliance. To maintain compliance with FCC RF exposure compliance requirements, please follow operation instructions as documented in this manual. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. The availability of some specific channels and/or operational frequency bands are country dependent and are firmware programmed at the factory to match the intended destination. The firmware setting is not accessible by the end user.

## I INSTRUCTION FOR REPORTING ADVERSE EVENT

MedWatch is the Food and Drug Administration's (FDA) program for reporting serious reactions, product quality problems, therapeutic inequivalence/failure, and product use errors with human medical products, including drugs, biologic products, medical devices, dietary supplements, infant formula, and cosmetics.

If you think you or someone in your family has experienced a serious reaction to a medical product, you are encouraged to take the reporting form to your doctor. Your health care provider can provide clinical information based on your medical record that can help FDA evaluate your report.

However, we understand that for a variety of reasons, you may not wish to have the form filled out by your health care provider, or your health care provider may choose not to complete the form. Your health care provider is NOT required to report to the FDA. In these situations, you may complete the Online Reporting Form yourself.

You will receive an acknowledgement from FDA when your report is received. Reports are reviewed by FDA staff. You will be personally contacted only if we need additional information.

## I SUBMITTING ADVERSE EVENT REPORTS TO FDA

1. Report Online at [www.accessdata.fda.gov/scripts/medwatch/index.cfm?action=reporting.home](http://www.accessdata.fda.gov/scripts/medwatch/index.cfm?action=reporting.home)
2. Consumer Reporting Form FDA 3500B. Follow the instructions on the form to either fax or mail it in for submission. For help filling out the form, see MedWatchLearn. The form is available at [www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm349464.pdf](http://www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm349464.pdf)
3. Call FDA at 1-800-FDA-1088 to report by telephone Reporting Form FDA 3500 commonly used by health professionals. The form is available at [www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm163919.pdf](http://www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm163919.pdf)

# WARRANTY INFORMATION

## TERMS & POLICY

**Nekteck Corporation ("Nekteck") warrants this product to the original purchaser to be free from defects in material and workmanship, under normal use and conditions, for a period of 1 year from the date of original purchase. Nekteck agrees, at our option during the warranty period, to repair any defect in material or workmanship or furnish an equal product in exchange without charge, subject to verification of the defect or malfunction and proof of the date of purchase. There is no other express warranty. This warranty does not apply:**

If the product has been modified from its original condition;

- If the product has not been used in accordance with directions and instructions in the user manual;
- To damages or defects caused by accident, abuse, misuse or improper or inadequate maintenance;
- To damages or defects caused by service or repair of the product performed by an unauthorized service provider or by anyone other than Nekteck;
- To damages or defects occurred during commercial use, rental use, or any use for which the product is not intended;
- To damages or defects exceeding the cost of the product.

1. Nekteck will not be liable for indirect, incidental, or consequential damages in connection with the use of the product covered by this warranty.
2. This warranty extends only to the original consumer purchaser of the product and is not transferable to any subsequent owner of the product regardless of whether the product is transferred during the specified term of the warranty.
3. This warranty does not extend to products purchased from unauthorized sellers. Nekteck's warranty extends only to products purchased from authorized sellers that are subject to Nekteck's quality controls and have agreed to follow its quality controls.
4. All implied warranties are limited to the period of this limited warranty.
5. This warranty gives you specific legal rights and you may also have other rights which vary from state to state. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

6. If you discover that your product is defective within the specified warranty period, please contact Customer Support via [support@Nekteck.com](mailto:support@Nekteck.com). DO NOT dispose of your product before contacting us. Once our Customer Support Team has approved your request, please return the product with a copy of the invoice and order ID.
7. Every Nekteck product automatically includes a 1-year warranty. To make the customer support process quick and easy, register your product online at [www.Nekteck.com](http://www.Nekteck.com).

**CORRECT DISPOSAL OF THIS PRODUCT (WASTE ELECTRICAL & ELECTRONIC EQUIPMENT)**

1. This marking shown on the product indicates that it should not be disposed with other house hold waste at the end of its life.
2. To prevent potential harm to the environment or to human health, please separate this product from other types of wastes and recycle it responsibly.
3. When disposing this type of product, contact the retailer where product was purchased or contact your local government office for details regarding how this item can be disposed in an environmentally safe recycling center.
4. Business users should contact their supplier and check the terms and conditions of the purchasing agreement. This product should not be mixed with other commercial wastes for disposal. This product is free of hazardous materials.

## I CUSTOMER SERVICE

1. If you have any product questions or concerns, please don't hesitate to contact us directly at [support@nekteck.com](mailto:support@nekteck.com) or visit our website at: [support@nekteck.com](http://support@nekteck.com).
2. We have a hassle-free warranty and customer service team for satisfactory solutions to any issue you might have.



**Email: [support@nekteck.com](mailto:support@nekteck.com)**

**We will reply to your message within 24 hours!**