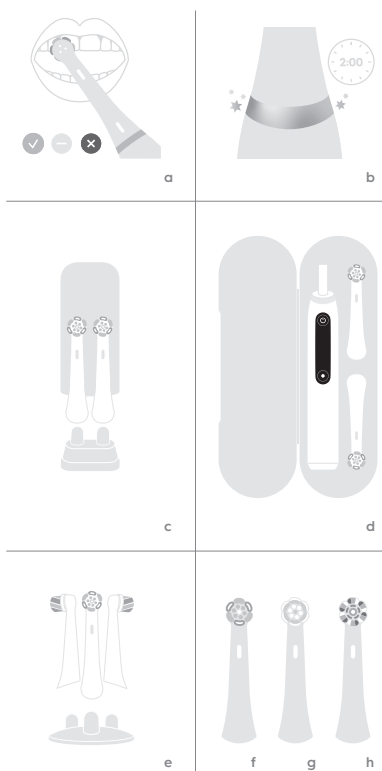
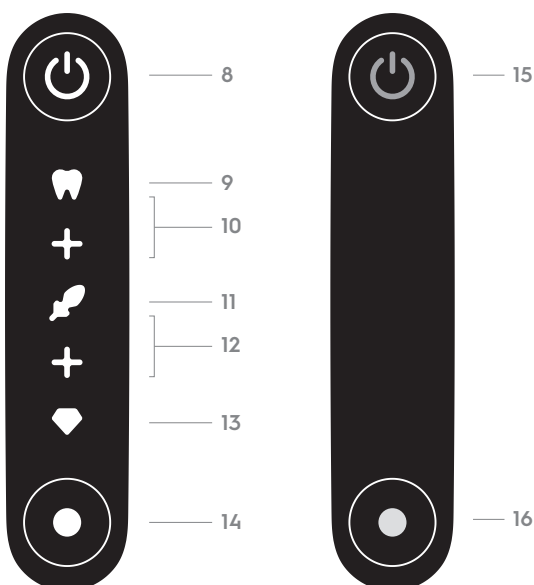
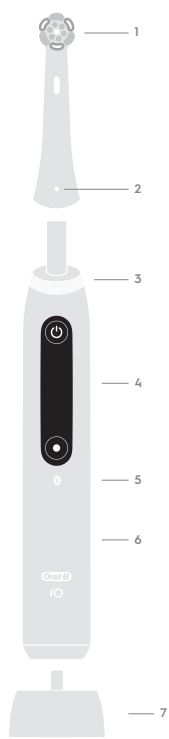
English

Description

- 1 Brush head
- 2 Brush head identification symbol
- 3 Smart Pressure Sensor / Light Ring
- 4 Brushing mode display
- 5 *Bluetooth*®
- 6 Handle
- 7 Charger
- 8 On/off button
- 9 «Daily Clean» mode
- 10 «Intense» mode
- 11 «Sensitive» mode
- 12 «Super Sensitive» mode
- 13 «Whiten» mode
- 14 Mode switch button
- 15 Low charge indicator
- 16 Brush Replacement Reminder

Features & Accessories (depending on model):

- a Smart Pressure Sensor
- b 2-minutes brushing reward illumination
- c Brush head holder with cap
- d Travel case
- e Brush head holder stand
- f Brush head «Ultimate Clean»
- g Brush head «Gentle Care»
- h Brush head «Radiant White»

Oral-B App:

- i Download Oral-B App and connect to *Bluetooth* (see back page)

Note: Content may vary based on model purchased.

Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

IMPORTANT

1. Periodically check the entire product/cord/accessories for damage. A damaged or non-functioning unit should no longer be used. If the product/accessories is damaged, take it to an Oral-B Service Center. The supply cord cannot be replaced. If the cord is damaged the appliance should be discarded. Do not modify or repair the product. This may cause fire, electric shock or injury.
2. Never insert any object into any opening of the appliance.
3. Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
4. For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
5. Cleaning and user maintenance shall not be made by children without supervision.
6. Children shall not play with the appliance.
7. Use this product only for its intended use as described in this manual. Do not use attachments (e.g. brush head refills) which are not recommended by the manufacturer.
8. Only use charger provided with the appliance.

WARNING

1. Do not place the charger/charging stand, power supply or charging travel case in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
2. This appliance contains batteries that are non-replaceable. The appliance has no user serviceable parts. Do not open and disassemble device except when disposing of the battery. Follow handling and disposal instructions of the battery/device described in this manual.
3. When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
4. This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution. For hygienic reasons, use individual brush handles and brush heads for each family member.
5. Small parts may come off, keep out of reach of children.
6. To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:
 - Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
 - If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
 - Replace brush head every 3 months or sooner if brush head becomes worn.
 - Clean brush head properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

Important Information

- This toothbrush and pouch (accessory) contains NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your Oral-B iO toothbrush.
- Your toothbrush can be used with your smartphone (find details under «Connecting your toothbrush to your Smartphone»). To avoid electromagnetic interference and/or compatibility conflicts, disable *Bluetooth* on your toothbrush handle (6) before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals.
- Persons with pacemakers should always keep the toothbrush more than 15 centimeters (six inches) from the pacemaker while turned on. Anytime you suspect interference is taking place disable *Bluetooth*.

Specifications

For voltage specifications refer to the bottom of the charging unit.

Using your toothbrush

Charging and Operating

- To charge the battery place the toothbrush on the plugged-in charger (7).
- During charging the on/off button is blinking white. Once it is fully charged the indicator light turns off. A full charge typically takes 16 hours.
- If the battery is running low the low charge indicator (15) blinks in red. Once the battery is empty, the motor stops. To charge the battery please put the brush on the charger.
- The toothbrush can be stored on the plugged-in charger to maintain it at full charge; a battery overcharge is prevented by the device. However, for environmental purposes we recommend to keep

charger in the power socket during charging only.

- The toothbrush switches off when exceeding the outside operating limits and/or standard battery current. To enable operation the toothbrush needs to be put on the charger.

Caution: This toothbrush contains a Li-Ion battery. Do not expose toothbrush to temperatures higher than 50°C / 122°F.

Note: Store the toothbrush at room temperature for optimal battery maintenance.

Start Brushing

Place the brush head (1) on the handle (6). A small gap should remain between handle and brush head. Activate the brushing mode display (4) by picking-up the brush, lifting it off the charger or pressing the on/off button (8).

Note: The brushing mode display turns off after some time. You may wake up the toothbrush to activate the display again.

Wet brush head and apply toothpaste. To avoid splashing, put the brush head into your mouth before switching the toothbrush on by pressing the on/off button.

Guide the brush head from tooth to tooth. Brush the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head (g)). Availability of modes may vary depending on model purchased.

Timer

The timer signals with a short stuttering sound in 30-seconds intervals to move to the next quadrant of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

The timer memorizes the elapsed brushing time, even when the toothbrush is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the toothbrush on the plugged-in charger.

Brushing reward illumination

The completion of a 2-minutes brushing event will be rewarded with a colorful light ring illumination (b).

Note: You may also customize the timer settings via Oral-B App.

Smart Pressure Sensor/ Light Ring (a)

Your toothbrush is equipped with a Smart Pressure Sensor (3) that has various functionalities indicated via different colors of light (a):

White = Low pressure indicator / default setting

Green = optimal pressure indicator

Red = high pressure indicator

Brushing with the recommended pressure is indicated via a green light.

If too much pressure is applied, a red light reminds you to reduce pressure. If the Smart Pressure Sensor illuminates in white (or your pre-selected color), we recommend applying more pressure.

When high pressure is applied the brush head movement will be reduced (in «Daily Clean», «Intense» and «Sensitive» modes).

Periodically check the operation of the Smart Pressure Sensor by pressing moderately on the brush head during use.

Note: You can replace the white default light ring color with a choice of individual colors via the Oral-B App.

Brushing modes (depending on model)

Your toothbrush offers different brushing modes indicated on the display:

«Daily Clean» (9) – Standard mode for daily cleaning

«Intense» (10) – For an extraordinary clean feeling (operated with a higher speed)

«Sensitive» (11) – Gentle, yet thorough cleaning for sensitive areas

«Super Sensitive» (12) – Super gentle cleaning for extra sensitive areas

«Whiten» (13) – Polishing for occasional or everyday use

Note: The sequence and availability of brushing modes may vary depending on model purchased.

Mode setting:

Wake up your iO toothbrush to identify selected mode for your brushing session. To change to others modes, successively press the mode button (14) during use. The toothbrush automatically starts in the mode selected last.

Note: You may also customize mode settings via Oral-B App.

Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads (f/g/h) that only fit your Oral-B iO toothbrush handle. These unique brush heads are dentist-inspired to leave your mouth feeling clean and refreshed.

When to exchange brush head

Some iO Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the bristle color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may be exerting too much pressure on teeth and gums. In general, we recommend replacing the brush head every 3 months.

Brush Replacement Reminder

The toothbrush is equipped with a «Brush Replacement Reminder» indicating the need to exchange the brush head approximately after 3 months when used twice per day. A yellow* light - accompanied by a vibration signal - will illuminate in the mode button (16).

To reset the «Brush Replacement Reminder» hold the mode button pressed for 2 seconds. A vibration and illumination signal confirms the reset. If no action is taken, an automatic reset will happen after several days.

*may appear orange

Note: You may also customize the brush replacement settings via the Oral-B App.

Connecting your toothbrush to your Smartphone/ Oral-B App

The Oral-B App offers various features:

- **Personalize your toothbrush settings** (e.g. choose your preferred light ring color, adjust stutter timer, mode sequence or brush replacement reminder settings)

- **Track your brushing behavior** (e.g. brushing time and applied pressure)
- **Real time brushing coaching** to help you brush every area for the right amount of time.

The Oral-B App is available for smartphones running iOS or Android™ operating systems. It can be downloaded for free from App Store^(SM)/Google Play™ or you may scan the QR code on the back page of this manual (i).

- Launch the Oral-B App. It will guide you through the entire *Bluetooth* pairing procedure. You may also find it under the App's settings (More -> Your brush).

Note: Oral-B App functionality is limited if *Bluetooth* is deactivated on your smartphone (for directions please refer to your smartphone's user manual).

- First pairing is indicated via the blue Smart Pressure Sensor light. Any App instructions will be displayed on your smartphone.
- First pairing needs to be done via Oral-B App. The toothbrush is automatically connected with your smartphone for future brush events.
- *Bluetooth* is enabled on your toothbrush by default. You may enable/ disable *Bluetooth* at any time by pressing the on/off button for 5 sec. Activation of *Bluetooth* is indicated via the light ring by a blue fade-in illumination. Deactivation is indicated by a blue fade-out illumination.
- Keep your smartphone in close proximity during brushing event. Make sure the smartphone is safely positioned in a dry place.


Note: Your smartphone must support *Bluetooth* 4.2 (or higher) / *Bluetooth* Smart in order to pair your toothbrush. Check app.oralb.com for compatible smartphones and country availability of the Oral-B App.

Caution: Follow your smartphone's user manual to ensure your phone/charger is designed for use in a bathroom environment.

Cleaning Guidelines

After each brushing, rinse your brush head under running water with the handle switched on. Remove the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle while the handle is switched on, as well as inside the brush head. Wipe all parts dry and let them air-dry before reassembling your toothbrush again.

Unplug the charger before cleaning. Clean the charger regularly with a damp cloth. Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.

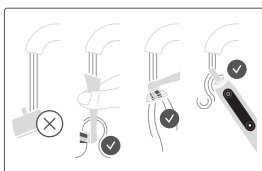Accessories (depending on model)

Travel Case (d)

Store the clean and dry toothbrush/brush heads in the travel case (d). The outer and inner case should be cleaned with a damp cloth only.

Note: Clean stains on travel case surface immediately. Store travel case in a clean and dry place.

Brush head holder stand (e) / Brush head holder with cap (c)

You may store your cleaned brush heads on the pins of the brush head holders. Make sure brush heads are dried before placing them in the holder and/or closing the cap. Brush head holders and cap are dishwasher-safe.

Zipper Case / Magnetic Pouch

Store the clean and dry toothbrush/brush heads in the Zipper Case / Magnetic pouch.

The outer and inner parts can be cleaned under running water. Store the clean and dry toothbrush along with the brush head holder (c) in the magnetic pouch. Make sure brush heads are dried before placing them into the holder and/or closing the cap. The outer and inner surfaces can be cleaned under running water. Store pouches in a clean and dry place.

Note: Clean stains on pouch surfaces immediately.

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your local region.

Battery removal at the end of the appliance's useful life

Do not disassemble the appliance except when disposing of the battery. Run down battery before removing. When taking out the battery for disposal, use protective clothing (gloves, safety glasses, etc.) to protect skin from a possible leak. Use caution not to short the positive (+) and negative (-) terminals. Insulate terminals (e.g. with tape) and put the battery in a plastic bag for disposal.

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to

provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes. Please also read our additional warranty statements: Warranty replacement brush heads.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts and accessories: Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative:

Please call 1-800-566-7252.

FCC Supplier's Declaration of Conformity

Responsible Party:	The Procter & Gamble Co. 1 P&G Plaza Cincinnati, OH 45202
Contact information:	www.oralb.com
Toothbrush:	Type 3776

FCC Compliance Statement: This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Section 15.21 Information to user: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Section 15.105 (b) Note: This equipment has been tested and found to comply with the Limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

Please also read our additional warranty statements: Warranty replacement brush heads on page 7.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B.

To obtain original Oral-B replacement parts & accessories:

Contact an authorized Oral-B Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To speak to an Oral-B Service Representative:

Please call 1-800-566-7252.

IC Compliance Notice: This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Satisfaction Guaranteed or Your Money Back!

Try it with no risk for 60 days!

If you are not 100% satisfied, call 1-877-593-4199 to receive a full refund. Return the product in its original packaging along with the original cash receipt postmarked within 60 days of purchase to receive a prepaid card in the amount of purchase. Terms of prepaid card apply. Your Oral-B toothbrush must be returned in its original packaging and be postmarked no later than 60 days after the date of purchase. Please allow 6–8 weeks for handling. In the US, the shipment cost to return the product is the responsibility of the consumer.

Sorry: No refunds will be made for any Oral-B product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service.

Warranty – Exclusions and Limitations

This limited warranty applies only to the new products manufactured by or for The Procter & Gamble Company, its Affiliates or Subsidiaries («P&G») that can be identified by the Braun / Oral-B trademark, trade name, or logo affixed thereto. This limited warranty does not apply to any non-P&G product including hardware and software. P&G is not liable for any damage to or loss of any programs, data, or other information stored on any media contained within the product, or any non-P&G product or part not covered by this limited warranty. Recovery or reinstallation of programs, data or other information is not covered under this limited warranty.

This limited warranty does not apply to (i) damage caused by accident, abuse, misuse, neglect, misapplication, or non-P&G product; (ii) damage caused by service performed by anyone other than Braun or authorized Braun service location; (iii) a product or part that has been modified without the written permission of P&G.

Information

Bluetooth® Indicator

Although all indicated functions on the *Bluetooth* device are supported, Oral-B does not ensure 100% reliability in connection, and feature operation consistency.

Operation performance and connection reliability are direct consequences of each individual *Bluetooth* device, software version, as well as the operating system of those *Bluetooth* devices and company security regulations implemented on the device.

Oral-B adheres to and strictly implements the *Bluetooth* standard by which *Bluetooth* devices can communicate and function in Oral-B toothbrushes.

However, if device manufacturers fail to implement this standard, *Bluetooth* compatibility and features will be further compromised and the user may experience functional and feature-related issues. Please note that the software on the *Bluetooth* device can significantly affect compatibility and operation.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable toothbrush damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable toothbrush at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high-quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

Trouble Shooting

Problem	Possible reason	Remedy
APP		
Handle does not operate (properly) with Oral-B App.	1. Oral-B App is closed. 2. <i>Bluetooth</i> is not activated on smartphone. 3. <i>Bluetooth</i> is disabled on handle. 4. <i>Bluetooth</i> connection to smartphone has been lost. 5. Your smartphone does not support <i>Bluetooth</i> 4.2 (or higher)/ <i>Bluetooth</i> Smart. 6. Outdated Oral-B App.	1. Launch Oral-B App. 2. Activate <i>Bluetooth</i> on smartphone (described in its user manual). 3. Enable <i>Bluetooth</i> by pressing on/off button for 5 seconds. 4. Pair handle and smartphone again via App settings. Keep your smartphone in close proximity when using it with handle. 5. Your smartphone must support <i>Bluetooth</i> 4.2 (or higher) / <i>Bluetooth</i> Smart in order to pair handle. Check app.oralb.com for compatible smartphones and country availability of the Oral-B App. 6. Download latest Oral-B App.
Reset to factory settings.	Original functionality wanted.	Use App to perform factory reset

TOOTHBRUSH		
Toothbrush is not operating (during first use).	1. Battery has a very low charge level. 2. Operating outside standard temperature.	1. Charge at least for 30 min. 2. Always operate in standard operating temperatures (> 0°C and < 40°C).
Toothbrush operates a short stuttering sound after 2 minutes only or does not operate any stuttering sound anymore.	Timer has been modified/ deactivated via App.	Use the Oral-B App to change timer settings or perform factory reset (see App Trouble Shootings).
Toothbrush does not charge.	1. Battery is at a very low charge level; it may take up to 30 min until brushing mode display turns on. 2. Ambient temperature for charging is out of valid range (< 0°C and > 40°C). 3. Charger may not be (properly) connected.	1. Charge at least for 30 min. 2. Recommended ambient temperature for charging is 5°C to 35°C. 3. Ensure charger is connected to mains voltage.
Rattling/loud noise during brushing or turning on.	1. Brush head may not be placed correctly onto the brush handle. 2. Brush head may be worn.	1. Check proper fit of brush head on handle. Small gap remains in between. 2. Replace brush head.
Mode button lights up in yellow color.	Brush Replacement Reminder indicates that your brush might need to be exchanged.	Exchange Brush head and reset timer by holding mode button for 2 seconds.

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Braun GmbH