

86.00 mm

120.00 mm

Numiere

FEATHER- Glow Meets Comfort

MANUAL



COMPANY INFORMATION

At Numiere, we believe skincare is a journey, not about fixing flaws. Lines, spots, and texture are part of your story, and we celebrate them. With years of expertise in the skincare and beauty industries, we're redefining who has access to the best tools and technologies to support their unique journey.

Our mission is to make high-quality skincare accessible to everyone. Numiere's tools are simple, effective, and designed for use anywhere, by anyone. We're here to empower you with restorative rituals and advanced innovations that enhance your skin's natural renewal process—wherever you are in your journey.

GENERAL DEVICE INFORMATION


The Numiere Feather LED Eye Mask is an FDA-cleared medical device that uses light therapy to effectively treat wrinkles around the eyes (periobital area).

DESCRIPTION OF THE DEVICE

The Numiere Feather LED Eye Mask combines three powerful light therapy treatments in one innovative device:

- Red Light Therapy – Helps reduce fine lines and wrinkles while promoting collagen production.
- Near-Infrared Light Therapy –Penetrates deeply to improve skin texture and elasticity.
- Far-Red Light Therapy – Enhances overall skin rejuvenation and circulation for a refreshed look.

PRODUCT STRUCTURE

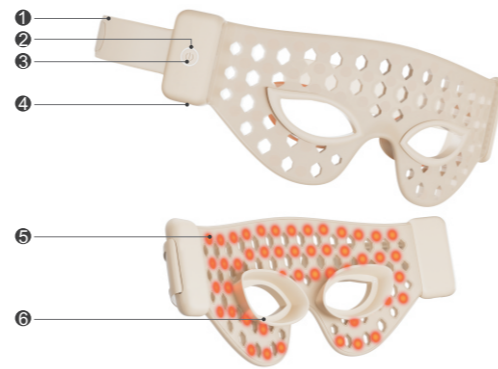

- 1 Bandage
- 2 Indicator light
- 3 Power/control key
- 4 Charging port
- 5 Light zone
- 6 Detachable Protective Eyewear


HOW TO USE :

- **Power On/Off:** Press and hold the power button for 1.5 seconds to turn the device on or off. Short press the button to adjust the gear.
- **Indicator Light:**When powered on, the indicator light will turn white, and the red light and NIR will activate simultaneously.The device will shut down in 10 minutes when there is no operation.
- **Low Battery Warning:**If the battery is low, the indicator light will flash quickly, and the device will automatically shut down after 5 seconds.
- **Charging Status:** During charging, the indicator light flashes slowly. Once fully charged, the indicator light will remain steady.

DEVICE SPECIFICATIONS

Energy:	max 70mw/cm2
LEDs:	Red:630nm +/-10nm
	Deep red:660nm+/-10nm
	NIR (Near-Infrared):850nm +/-10nm
LED quantity:	51 pcs
Program Time:	10 minutes
Power:	Max2.5w
Battery Capacity:	500mAh
Operating Conditions:	Operating Temperature:0°C to 40°C
Storage and Transport Temperature:	-5°C to 55°without condensation
Operating Humidity:	30% to 95%
Storage and Transport Humidity:	20% to 85% without condensation
Charging supply:	5V = 1A
Charging cable :	TYPE C
Charging cable length :	80cm
Charging time :	about 90 minutes
Time of use :	80minutes
Recommended use :	1 time per day,3-7times per week
Materials :	Silicone_PC,Nylon,PU

DEVICE ACCESSORIES

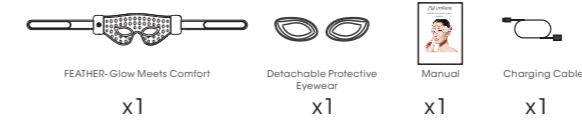

ADVERSE REACTION

If you have any adverse reaction from using this device,such as discomfort, skin discoloration, skin irritation, or other concerns, stop using it immediately and consult your health professional. Skin irritation, bruising, contact burns,dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your health professional.If you experience a persistent headache and/or otherpersistent painful sensations while using this device, stopusing it and consult with your health professional

CONTRAINDICATIONS

- DO NOT use the device if you have any active cancer or cancerous lesion.
- DO NOT use the device if you are pregnant, planning to become pregnant, or breast feeding.The device has NOT been tested on pregnant or breast-feeding women and therefore the risk to the fetus, neonate or pregnant women is unknown.
- DO NOT use the device if you are under the age of 22.
- DO NOT use the device over areas with phlebitis,thrombo phlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to silicone or plastics.The device contains plastics and silicones that come into contact with the skin during use.When in doubt of allergies, contact your health professional before using the device.
- DO NOT use the device if you suffer from any photosensitive disorder (sensitization to light) or
- DO NOT use the device if you suffer from any serious skin condition such as severe acne.
- DO NOT use the device if you are taking any medication or drugs that can cause photosensitivity. If you are not sure about your medication, contact your health professional.
- DO NOT use the device if you suffer from light induced headaches.
- DO NOT use the device if you have a medical history of seizures or epilepsy triggered by light.

- DO NOT use the device if you have the Fitzpatrick Skin in Type V and Type VI.
- DO NOT continuously press the ON/OFF button.
- Note:** We highly suggest being cautious when taking any (skin-related) medications. If in doubt, consult your health professional.
- Note:** If you are under the care of a physician, consult them before using this device.

WARNINGS

- **STOP** using the device and consult your health professional if you experience an adverse reaction from using the device.
- **STOP** using the device if you experience any discomfort or difficulty breathing or swallowing occurring from application of the device on the neck.
- **STOP** using the device if you experience symptoms of hypotension (dizziness or light-headedness) while using the device.
- **STOP** using the device if its overheating. Turn it off immediately.
- **STOP** your treatment if it feels hot or your skin feels uncomfortably warm.
- **DO NOT** let children play with the device. Cleaning and charging the device shall not be done by children.
- **DO NOT** use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- **DO NOT** use the device on spots without pigmentation or areas on the skin that have vitiligo.
- **DO NOT** shine the light from the LEDs directly into your eyes or someone else's eyes. This may cause temporary loss of sight.
- **DO NOT** place the device on any open wounds.
- **DO NOT** use the device outdoors to prevent damage to the product.
- **DO NOT** sit, lay down on, or step on the device. This will break the device.
- **DO NOT** drop or insert any object into any openings on the device.
- **DO NOT** use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs, contact customer support.

- **DO NOT** use or store the device near hot surfaces.
- **DO NOT** charge the device when it is close to or might come into contact with water.
- **DO NOT** submerge the device in water.
- **DO NOT** store the device in extreme heat or cold.
- **DO NOT** place the device in direct sunlight.
- **DO NOT** attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. It contains batteries that are non-replaceable.
- **DO NOT** perform maintenance on or take apart the product. In case of malfunctioning or damage, contact customer support.
- **DO NOT** use the device if it is damaged or if the charging cable is damaged.
- **DO NOT** use the device with non-manufacturer-approved accessories.
- **DO NOT** use solvents or abrasive cleaners on the device.
- **DO NOT** use the device in any other way than specified in this manual.
- **DO NOT** use the device if you have had a recent surgery or other surgical procedures.
- **DO NOT** use the device over a bacterial infection, as this may result in the worsening of symptoms.
- **DO NOT** use the device on painful areas. Consult your health professional when you encounter painful areas during use.
- **DO NOT** use the device on or near children. Do not let it be used by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned, or stored near children or those with reduced physical or mental abilities.
- **DO NOT** use the device with any USB cable other than the one provided by Numiere.
- **DO NOT** use the device over areas of skin that lack normalsensation.
- **DO NOT** use the device on normal, intact, clean, and healthy skin.
- **DO NOT** use the device while driving, operating machinery, or during an activity in which voluntary muscle contractions may put users or the public at any risk of injury.

- **DO NOT** use the device on swollen, infected, or inflamed areas or skin eruptions.
- **DO NOT** use the device in high humidity locations such as the shower or sauna.
- **DO NOT** exceed the recommended treatment time when using the device.
- If there has been sunburn on your skin in the past two weeks, it is recommended not to use this product.
- This device is intended to be used by one person, we recommend to avoid sharing with others.
- Note:** The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cable is damaged, discontinue use. Contact Numiere customer support to arrange a cable replacement.

PRECAUTIONS

- **DO NOT** scratch or damage the light-emitting surface.
- Not operating the device in accordance with the instructions in the manual will void your warranty.
- Read all instructions carefully before using this device.
- Store the device according to the protocol outlined in this manual.
- This device is intended to be used by one person only.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your health professional if you are unsure about your health status.
- The device's battery is not removable. If you are disposing of the device, do not burn it or throw it in the trash. Dispose of the device according to local regulations.
- Inspect the device and do not use it if you see any signs of damage, such as internal battery leaks. If anything seems out of the ordinary, contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cord is damaged, discontinue use. Contact customer support to arrange a replacement.
- Consult your health professional before using the device

- if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anesthetic, or desensitized areas.
- We highly suggest using caution if you are taking any medications. If in doubt, consult your health professional. If using the device causes persistent pain or discomfort, discontinue use immediately.
- There may be a perception of flashing lights during the use of the device, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your device, consult your health professional.

DISCLAIMER

This device is not intended to diagnose, cure, or prevent disease. Numiere makes no claims, representations, or warranties regarding the ability of this device to cure any physical, skin, or mental condition. A qualified health professional should be consulted regarding any condition requiring medical attention.

RISKS AND BENEFITS

POTENTIAL RISKS

- Immediately after treatment, skin may appear red. Skin should return to normal within 1-2 hours.
- Temporary headaches may develop immediately following the light treatment. Headaches can be treated with an appropriate pain killer and rehydration.

POTENTIAL BENEFITS

- Red Light Therapy: reduces the appearance of fine lines and wrinkles.
- Far-Red Light Therapy: improves skin wellness and deeply rejuvenates the skin.
- Near-Infrared Light Therapy: helps manage suboptimal local blood flow and circulation while addressing muscle and joint pain

POTENTIAL BENEFITS

We highly recommend users perform a light sensitivity test before their first use of the device.

1. Place the device on the inner forearm.
2. Turn on the device to start a treatment.
3. Wait until the treatment is done. The full treatment will take 3 minutes to complete.
4. If your skin feels hot or uncomfortable during the treatment, stop immediately! You might be sensitive to light.
5. The device will automatically turn off after the treatment is done. Remove the device and wait for at least 24 hours before your next treatment. Check your skin regularly over those 24 hours. If your skin appears red and/or itchy, the device is unsuitable for you. In case of sensitivity to the light, please contact customer support.

MANUFACTURE INFORMATION

Heat in a click llc
 Numiere
 1975 Tigertail blvd Dania FL 33004
 PH: (888) 691-1905
 www.numiere.com