

Thank you for your purchase!

Scan the code below
for help getting started



Got questions?

We're here to help!

For Customer Support, Please Visit
<https://mscustomercare.thethreesixtygroup.com>
or call

1-800-374-2744.

Available Monday – Friday, 9am – 5pm PST
2-Year Warranty


QUICK PROGRAM GUIDE

PULSE


Press to start default program.

Default


- Program is Pulse all 5 zones.
- Pressure level is 1.
- Run time is 10 min.

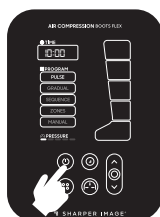


Press to add 10 min.
Max30 min.

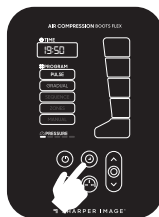

Press to raise pressure level from 1 to 5.

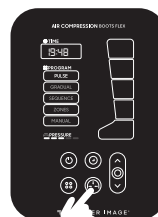


PULSE Program in action

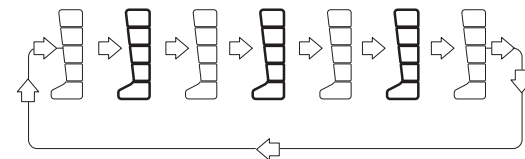


GRADUAL


Press to start default program.

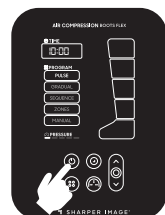


Press once to select the GRADUAL program.

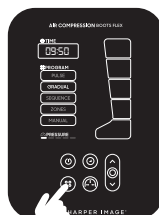


Press to add 10 min.
Max30 min.

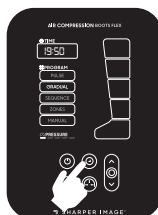

Press to raise pressure level from 1 to 5.

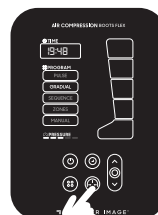


GRADUAL program in action

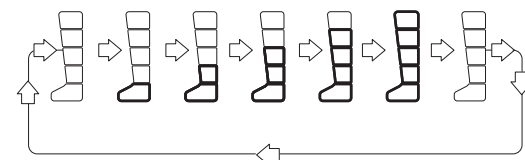


SEQUENCE


Press to start default program.

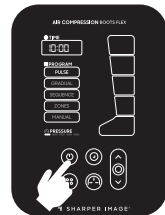


Press twice to select the SEQUENCE program.

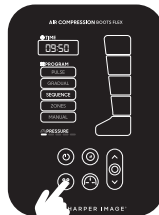


Press to add 10 min.
Max30 min.

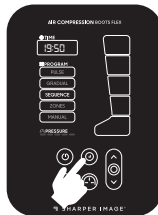

Press to raise pressure level from 1 to 5.

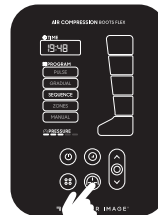


SEQUENCE program in action

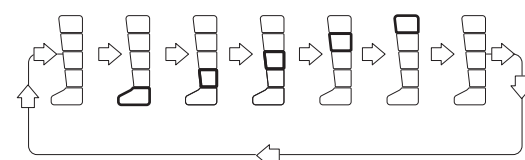


ZONES


Press to start default program.

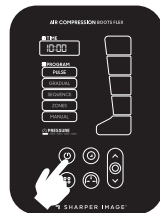


Press three times to select the ZONES program.

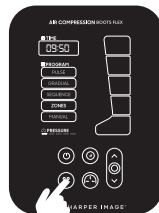


Press or arrows to choose a zone to include in the zones massage.

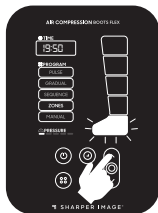

Press to confirm the zone. Select the next zone or and press to confirm.

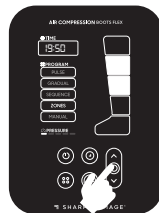


Long press to start the pulse pressure massage for the selected zones.

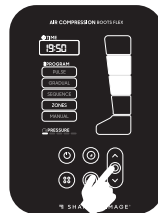

Press to add 10 min.


Press to raise pressure level.

ZONES program in action

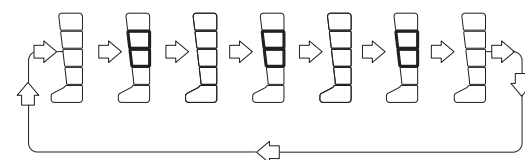


MANUAL


Press to start default program.

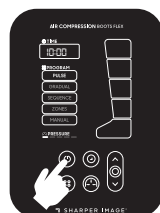


Press four times to select MANUAL program.

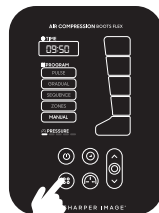

Press or arrows to choose a zone to include in the zones massage.

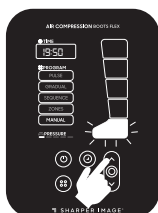

Press to confirm the zone. Select the next zone or and press to confirm.

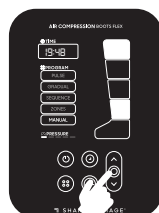

Long press to start the pulse pressure massage for the selected zones.



Press to add 10 min.

Press to raise pressure level.

MANUAL program in action

