

HIGHERDOSE

A close-up portrait of a woman with dark skin and short hair, looking slightly to the right. She is wearing a light blue neck enhancer device around her neck. The device has a glowing red and orange light around the neck area. The brand name 'HIGHERDOSE' is printed in white on the blue device. The background is dark and out of focus.

RED LIGHT NECK ENHANCER

INSTRUCTION MANUAL

TABLE *of* CONTENTS

WELCOME	03
THE HIGHERDOSE RED LIGHT NECK ENHANCER	04
MEDICAL WARNINGS	04
PRODUCT INFORMATION - SPECIFICATIONS	05
WHAT'S INSIDE THE BOX	06
SENSITIZATION TEST	06
POTENTIAL ADVERSE REACTION	06
CONTRAINDICATIONS & PRIOR MEDICAL ADVICE	07
WARNINGS	07
PRECAUTIONS	08
TRAVEL AND INTERNATIONAL USE	09
DISCLAIMER	09
POTENTIAL RISKS	09
HOW TO USE YOUR RED LIGHT NECK ENHANCER	10
ALTERNATIVE WAYS TO USE	11
TREATMENT PROTOCOL	11
CHARGING	12
FREQUENTLY ASKED QUESTIONS	12
CARE AND STORAGE	13
CLEANING	13
TROUBLESHOOTING	14
WARRANTY	15

LET'S GET STARTED

RED LIGHT NECK ENHANCER INSTRUCTION MANUAL



WELCOME to your NEW HIGHERDOSE RED LIGHT NECK ENHANCER

READY TO IGNITE
NATURALLY GLOWING SKIN?

Before using your **RED LIGHT NECK ENHANCER** read this **user manual thoroughly** to understand correct product usage, view important safety precautions, and learn how to boost the benefits of your session.

Save this manual for future reference.

The **HIGHERDOSE RED LIGHT NECK ENHANCER** has been carefully engineered for effectiveness and safety. Basic safety precautions should always be observed when using this product. Before using this product, please read this manual carefully.

HIGHERDOSE RED LIGHT NECK ENHANCER

The HIGHERDOSE RED LIGHT NECK ENHANCER is an at-home wearable LED phototherapy device with two proven wavelengths of light: 630 nm (Red Light) and 830 nm (near Infrared Light). Both of these lights are known to rejuvenate the skin.

This red light device has been designed to fit your neck and décolleté so that the maximum amount of light hits the targeted skin. We've included multiple straps to ensure a secure, comfortable fit and to offer a wide range of treatment area options. The HIGHERDOSE RED LIGHT NECK ENHANCER covers wider areas of the chest and around the neck offering more LED bulbs per device than others on the market.

MEDICAL WARNINGS

Consult your physician before using this product. Do not use The HIGHERDOSE RED LIGHT NECK ENHANCER if you are taking any medication that can cause photosensitivity.


Those with any skin sensitivities or any medical conditions should consult their doctor before use. In rare cases, prolonged exposure to heat may result in a condition causing discoloration of the skin known as Erythema Ab Igne ("Toasted Skin Syndrome"). If you notice any skin discoloration, discontinue use immediately.

PRODUCT INFORMATION - SPECIFICATIONS

LED WAVELENGTHS	RED: 630nm 26 mW / cm2 & NIR 830nm 24mW / cm2
TOTAL IRRADIANCE	50mW /cm2
TOTAL LED / DIODES	75 dual-core LEDs (150 DIODES)
OPERATING TEMPERATURE	0°C to 40°C 32°F to 104°F
STORING AND TRANSPORT TEMPERATURE	-5°C to 55°C 23°F to 131°F
OPERATING HUMIDITY	30% to 95% without condensation
STORAGE AND TRANSPORT HUMIDITY	20% to 95% without condensation
CHARGING SUPPLY	5V = 1A
CHARGING PORT	USB-C
CHARGING CABLE	USB-C to USB-A
CHARGING CABLE LENGTH	3.28 ft / 1.2 meter
BATTERY	3.7V, 2600 mAh
BATTERY TYPE	Lithium-ion
CHARGING TIME	3.5-4 hours
TREATMENT TIME	10 or 20 minutes
DIMENSIONS OF THE DEVICE	332 X 310 X 5mm
WEIGHT	14.8 oz / .422 kg
ELECTRICAL OUTPUT	5V/1A/2600mAh
ELECTRICAL OUTPUT FROM POWER ADAPTER	100-240V

WHAT'S INSIDE THE BOX?

1. RED LIGHT NECK ENHANCER
2. Controller
3. Power Supply and Country-Specific Adaptors
4. USB C to USB A connector
5. Adjustable Straps
6. User Manual (not pictured)
7. Goggles

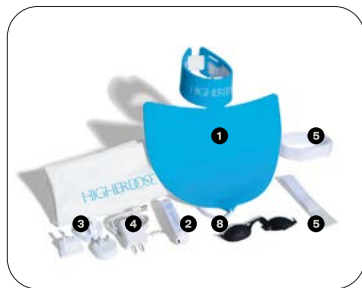

SENSITIZATION TEST

Before you begin the first treatment, we recommend the sensitization test on your skin. Complete the following steps to conduct the test:

- Pick one of the inner arms and place the device on the desired area.
- Press POWER and allow the device to run for a 10-minute treatment program. If your skin feels hot or uncomfortable STOP using the RED LIGHT NECK ENHANCER immediately.
- Initially, your skin might look a little red. Wait 3-4 hours after the treatment. If the skin on your inner arm still appears red or itchy, this product might not be suitable for you. Please contact HIGHERDOSE customer service team at help@higherdose.com
- A pair of goggles is provided. If you have photosensitive eyes, you should use them during the treatment.

POTENTIAL ADVERSE REACTION

If you have any adverse reaction, such as discomfort, skin discoloration, skin irritation, or other concerns, from using this device, stop using it immediately and consult your health professional.

Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your doctor.

CONTRAINDICATIONS & PRIOR MEDICAL ADVICE

If you have any questions regarding your current medical condition and the use of The HIGHERDOSE RED LIGHT NECK ENHANCER, please consult your medical doctor before use.

If any of the following conditions apply, seek permission from your doctor prior to use:

- Taking medications
- Under 12 years old
- Chronic Conditions

If any of the following conditions apply, The HIGHERDOSE RED LIGHT NECK ENHANCER's use is contraindicated:

- Lupus Erythematosus
- Photosensitive eczema
- Albinism
- Photosensitive Disorders (sensitivity to light)
- Sensitivity to Heat

This product does not diagnose, treat, or cure any disease, nor is it implied or inferred to cure or treat any disease. In the rare event that you experience pain and/or discomfort while using this product, discontinue use immediately and consult your doctor.

WARNINGS

- **STOP** using the device and consult your doctor if you experience an adverse reaction from using the device.
- **STOP** using the device if you have symptoms of hypotension (dizziness or light-headedness) while using the device.
- **STOP** using the device if it is overheating. Turn off the device immediately.
- **STOP** your treatment if it feels hot or your skin feels uncomfortably warm.
- **DO NOT** let the children play with the device. Cleaning and charging the device shall not be done by children.
- **DO NOT** use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- **DO NOT** use the device on spots without pigmentation or areas on the skin that have vitiligo.

- **DO NOT** use on any open wounds.
- **DO NOT** use the device outdoors to prevent damage to the product.
- **DO NOT** sit, lay down, or step on the device. This will break the device.
- **DO NOT** drop or insert any object into any openings of the device.
- **DO NOT** use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs contact customer support at help@higherdose.com.
- **DO NOT** use or store the device near hot surfaces.
- **DO NOT** charge the device when it is close to or might get in contact with water.
- **DO NOT** submerge the device in water.
- **DO NOT** place the device in direct sunlight.
- **DO NOT** attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. This device contains batteries that are non-replaceable.
- **DO NOT** perform maintenance or take apart the product. In case of malfunctioning or damage, contact customer support.
- **DO NOT** use the device if damaged or when the charging cable is damaged.
- **DO NOT** use the device with non-manufacturer-approved accessories.
- **DO NOT** use solvents or abrasive cleaners on the device.
- **DO NOT** use the device in any other way than specified in this manual.
- **DO NOT** use the device if you have had a recent surgery or other surgical procedures.
- **DO NOT** use the device over a bacterial infection as this may result in the worsening of symptoms.
- **DO NOT** use the device on painful areas. Consult your doctor when you experience painful areas during use.
- **DO NOT** use the device on, near, or by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned, or stored near children or those with reduced physical or mental abilities.
- **DO NOT** use the device with any other USB-C cable than provided by HIGHERDOSE.
- **DO NOT** use the device with any topical other than specified by HIGHERDOSE.

Note: The device should only be supplied with 5V DC low voltage. If any part of the main cable is damaged, discontinue use. Contact customer support at help@higherdose.com to arrange a cable replacement.

PRECAUTIONS

- Read all instructions carefully before using this device
- Do not scratch or damage the light-emitting surface.
- It will void your warranty if you do not operate the device in accordance with the instructions in the manual.
- Store the device according to the protocol in this manual.
- This device is intended to be used by a single person.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your doctor if you are unsure about your health status.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put it into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of damage such as internal battery leaks. Contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main cable is damaged, discontinue use. Contact customer support (help@higherdose.com) to arrange a cable placement.
- Consult your doctor before using the device if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anesthetic, or desensitized areas.
- If using any medications, consult your doctor before using the product.

TRAVEL AND INTERNATIONAL USE

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are traveling to a different region than where you bought the device, please ensure that the device is only supplied with a 5V DC low voltage.
- Confirm with the steward that the device can be carried and/or used on the airplane prior to flight.

DISCLAIMER


DO NOT use the device to treat any other conditions apart from the conditions that are indicated.

POTENTIAL RISKS

- Immediately after treatment, the skin may appear red. The skin should return to normal within 1-2 hours.
- Temporary headaches may appear after the light treatment. Consult your health professional if the headaches persist over a longer period.

HOW TO USE YOUR RED LIGHT NECK ENHANCER

1. Before starting, make sure the battery is fully charged and you are using the right size strap for the area you want to treat. See the images on pages 20-21 for treatment area alternatives.
2. Start with freshly cleansed skin. Make sure to remove all make-up, dirt, and daily pollutants from your neck, chest, and décolleté area.
3. If you have our LIGHT-ACTIVATED GLOW SERUM, now is the perfect time to apply 2-3 pumps to your neck and décolleté. Do not apply sunscreen before using this device.
4. Plug your handheld charger into the USB plug located at the bottom of the cord on your HIGHERDOSE RED LIGHT NECK ENHANCER.
5. Secure the strap around your neck so that it is comfortable and close to your skin.
6. Press and hold the ON/OFF button until the green light turns on. Press once for 10 minutes and twice for 20 minutes. The timing buttons will illuminate for the treatment time you've chosen. You may stop the treatment at any time by pressing the ON/OFF button.
7. Now it's time to soak up the red light! Relax and enjoy the rejuvenating benefits of red light therapy. The cordless design allows you to move around effortlessly.
8. We suggest reapplying a little LIGHT-ACTIVATED GLOW SERUM post-treatment or following up with your favorite skin care cream.
9. Clean the device using a water-based wipe or a damp clean cloth. Make sure the device is NOT connected to the power supply when cleaning.

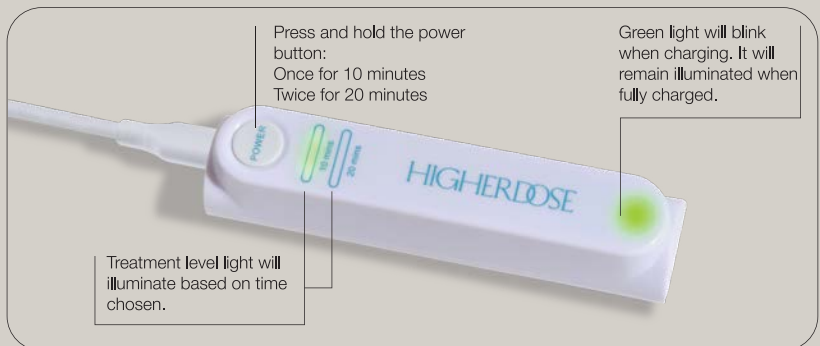

ALTERNATIVE WAYS TO USE YOUR RED LIGHT NECK ENHANCER

Target your neck and décolleté with more coverage area than competitors. Increased neck length allows for full neck coverage all around.



Target neck/back

Alternative straps allow for full chest area coverage



TREATMENT PROTOCOL

Use the device 4-5 times a week for a minimum of 4 weeks. After that, you can continue using your RED LIGHT NECK ENHANCER 2-4 times a week based on your lifestyle and skin care goals. Consistency is key and the results are cumulative, you will continue to see results over time with consistent use. It is not recommended to use the device more than 7 times per week.

CHARGING

Charge your HIGHERDOSE device fully for 3.5- 4 hours before your first use. Only use a 5V adapter or computer USB port (avoid fast chargers as they can damage the battery). Your battery will come partially charged.

To charge the device, insert the USB-C cable into the controller. Insert the USB-A side of the cable into a USB power supply. When the controller is charging, the LED indicator on the device flashes green, and when the controller is fully charged, the green light will remain on.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use. The device will not work while being charged.

FREQUENTLY ASKED QUESTIONS

WHAT IS RED LIGHT THERAPY?

- Red Light Therapy is a process where light-emitting diodes (LED) are used to deliver red (630 nm) and near-infrared light (830 nm) on your face or any other body parts. Red Light provides wavelengths of light similar to the sunlight but without the heat or damaging UV rays.

WHAT IS THE NECK ENHANCER MADE OF?

- The outer case of the RED LIGHT NECK ENHANCER is made with soft and flexible, medical-grade silicone. The LED bulbs are attached to this.

WHAT ARE THE BENEFITS?*

- Reduces redness- Low wavelength red light targets uneven skin tone.
- Boosts mood- Red light and near-infrared relax the body and mind
- Rejuvenates skin- Red light therapy helps to revitalize skin and enhance natural glow.

HOW MANY USES DO I NEED BEFORE I START SEEING BENEFITS?*

- It is possible to feel refreshed and invigorated after just one session, and many people see benefits after just one or two uses. The results may also continue to build over continued use.

***These statements have not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease.*

However, if you have photosensitivity we advise you to exercise caution and common sense. Listen to your body and do what is most comfortable.

CARE AND STORAGE

- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in dry conditions.
- Store the device in the supplied packaging. Make sure the device will not bend, stretch, or crush while being stored or when traveled with.
- Avoid storing the device in dusty conditions, in direct sunlight, or in damp and humid conditions.
- This device contains no serviceable parts.
- Before each use, with the device switched off, clean the device following the cleaning instructions.

IMPORTANT

- The product has a limited lifecycle of two years under nominal use. Following the treatment protocol for two years is considered nominal use.
- The device must be handled with care.

CLEANING

- Clean the device with a water-based wipe or a damp clean cloth. Try our HIGH MAINTENANCE CLEANER for a non-toxic, probiotic, and natural ingredient gentle cleaner. Safe for you and the environment.
- You can use a non-toxic cleaner, make-up removal wipe or baby wipe to clean the RED LIGHT NECK ENHANCER.

- DO NOT clean the device while it is turned on.
- DO NOT clean the device while charging.

TROUBLESHOOTING

If the actions suggested below fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please contact the customer service team at help@higherdose.com.

- DO NOT attempt to investigate any further.
- DO NOT attempt to modify the device in any way.
- The device has NO parts you can repair yourself

PROBLEM	CAUSE	SOLUTION
The device doesn't switch on	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the cable into the device correctly and charge the battery until the green light stays illuminated.
	The power outlet didn't work during charging.	Try a different power outlet.
	The on/off button has not been pushed for 1 full second.	Push the on/off button for 1 second straight.
The device switches off unexpectedly.	The battery level is too low.	Charge the device.

HIGHERDOSE RED LIGHT NECK ENHANCER

One (1) Year Limited Warranty Updated May 16th, 2024

IMPORTANT: BY USING YOUR HIGHERDOSE RED LIGHT NECK ENHANCER YOU ARE AGREEING TO BE BOUND BY THE TERMS OF THE HIGHERDOSE ONE (1) YEAR LIMITED WARRANTY ("WARRANTY") AS SET OUT BELOW.

DO NOT USE YOUR PRODUCT UNTIL YOU HAVE READ THE TERMS OF THE WARRANTY. IF YOU DO NOT AGREE TO THE TERMS OF THE WARRANTY, DO NOT USE THE PRODUCT AND RETURN IT WITHIN THE RETURN PERIOD STATED IN HIGHERDOSE'S RETURN POLICY FOR A REFUND.

HOW CONSUMER LAW RELATES TO THIS WARRANTY

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, BUT YOU MAY HAVE OTHER RIGHTS THAT VARY FROM STATE TO STATE (OR BY COUNTRY OR PROVINCE). FOR A FULL UNDERSTANDING OF YOUR RIGHTS YOU SHOULD CONSULT THE LAWS OF YOUR COUNTRY, PROVINCE, OR STATE.

WARRANTY LIMITATIONS SUBJECT TO CONSUMER LAW

TO THE EXTENT PERMITTED BY LAW, THIS WARRANTY AND THE REMEDIES SET FORTH ARE EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, REMEDIES, AND CONDITIONS, WHETHER ORAL, WRITTEN, STATUTORY, EXPRESS OR IMPLIED. HIGHERDOSE DISCLAIMS ALL STATUTORY AND IMPLIED WARRANTIES, INCLUDING WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS, TO THE EXTENT PERMITTED BY LAW. IN SO FAR AS SUCH WARRANTIES CANNOT BE DISCLAIMED, HIGHERDOSE LIMITS THE DURATION AND REMEDIES OF SUCH WARRANTIES TO THE DURATION OF THIS EXPRESS WARRANTY AND, AT HIGHERDOSE'S OPTION, THE REPAIR OR REPLACEMENT SERVICES DESCRIBED BELOW. SOME STATES (COUNTRIES AND PROVINCES) DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY (OR CONDITION) MAY LAST, SO THE LIMITATION DESCRIBED ABOVE MAY NOT APPLY TO YOU.

WHAT IS COVERED BY THIS WARRANTY?

HIGHERDOSE LLC of 42 Broadway, 12th floor, #212, New York City, New York, 10004 ("HIGHERDOSE") warrants The HIGHERDOSE- branded RED LIGHT NECK ENHANCER and The HIGHERDOSE- branded accessories contained in the original packaging ("HIGHERDOSE RED LIGHT NECK ENHANCER") against defects in materials and workmanship when used normally in accordance with HIGHERDOSE's published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser ("Warranty Period"). HIGHERDOSE's published guidelines include but are not limited to information contained in technical specifications, user manuals, and service communications.

WHAT IS NOT COVERED BY THIS WARRANTY?

This Warranty does not apply to any non-HIGHERDOSE branded products, even if packaged or sold with HIGHERDOSE-branded products. HIGHERDOSE does not warrant that the operation of The HIGHERDOSE RED LIGHT NECK ENHANCER will be uninterrupted or error-free. HIGHERDOSE is not responsible for damage arising from failure to follow instructions relating to The HIGHERDOSE RED LIGHT NECK ENHANCER use.

This Warranty does not apply:

- (a) To cosmetic damage, including but not limited to scratches, dents, and broken plastic on ports, cords, or controllers, unless failure has occurred due to a defect in materials or workmanship;
- (b) To damage caused by accident, abuse, misuse, fire, liquid contact, earthquake, or other external cause;
- (c) To damage caused by operating The HIGHERDOSE RED LIGHT NECK ENHANCER outside HIGHERDOSE's published guidelines;
- (d) To damage caused by service performed by anyone who is not a representative of HIGHERDOSE;
- (e) To a HIGHERDOSE RED LIGHT NECK ENHANCER that has been modified to alter functionality or capability without the written permission of HIGHERDOSE;

(f) To defects caused by authorized users of the product (eg. by presenting proof of purchase). normal wear and tear or otherwise due to the normal aging of The HIGHERDOSE RED LIGHT NECK ENHANCER;

(g) Damage caused by pets;

(h) If HIGHERDOSE receives information from relevant public authorities that the product has been stolen, and you cannot prove in any way that you are the authorized user of the product (eg, by presenting proof of purchase).

YOUR RESPONSIBILITIES

Before receiving warranty service, HIGHERDOSE or its agents may require that you furnish proof of purchase details, respond to questions designed to assist with diagnosing potential issues and follow HIGHERDOSE's procedures for obtaining warranty service.

Following warranty service your HIGHERDOSE RED LIGHT NECK ENHANCER or a replacement item will be returned to you, subject to applicable product model updates and availability.

Important: Do not attempt to repair The HIGHERDOSE RED LIGHT NECK ENHANCER, the Controller, or any other aspect of The HIGHERDOSE product. Attempting to repair The HIGHERDOSE product is dangerous and may cause damage that is not covered by this Warranty. Only HIGHERDOSE should perform service on the product.

NO COMMERCIAL USE

All HIGHERDOSE- branded consumer products including RED LIGHT NECK ENHANCER, are designed and intended only for home use. HIGHERDOSE- branded consumer products are not intended to be used in any commercial or institutional setting or application. These include but are not limited to commercial gyms, commercial spas, wellness centers, chiropractic offices, physician offices, any other medical facilities, churches, police/fire departments, government agencies, apartment/retirement/condominium communities, hotels, cruise ships, schools, or universities.

Using any HIGHERDOSE-branded product in a commercial or institutional setting will void the product's warranty and may open the owner/operator of the setting to liability.

WHAT WILL HIGHERDOSE DO IN THE EVENT THE WARRANTY IS BREACHED?

If during the Warranty Period, you submit a claim to HIGHERDOSE in accordance with this warranty, HIGHERDOSE will, at its option:

- Repair The HIGHERDOSE RED LIGHT NECK ENHANCER using new or previously used parts that are equivalent to new in performance and reliability.
- Replace the Controller.
- Replace The HIGHERDOSE RED LIGHT NECK ENHANCER with the same model (or with your consent a product that has similar functionality) formed from new and/or previously used parts that are equivalent to new in performance and reliability.
- Exchange The HIGHERDOSE RED LIGHT NECK ENHANCER for a new item.

HOW TO OBTAIN WARRANTY SERVICE?

You must contact HIGHERDOSE via email at help@HIGHERDOSE.com in order to seek warranty service.

A HIGHERDOSE representative will help determine whether your HIGHERDOSE RED LIGHT NECK ENHANCER requires service and, if it does, will inform you whether and how HIGHERDOSE will provide it.

WARRANTY SERVICE OPTIONS

HIGHERDOSE will provide warranty service through mail-in service. If HIGHERDOSE determines that your HIGHERDOSE RED LIGHT NECK ENHANCER is eligible for mail-in service, HIGHERDOSE will send you prepaid waybills and instructions on how to properly pack and address your HIGHERDOSE RED LIGHT NECK ENHANCER, so that you may ship your HIGHERDOSE RED LIGHT NECK ENHANCER to HIGHERDOSE. Instructions will be sent to you via email. Once service is complete, HIGHERDOSE will return the product to you. HIGHERDOSE will pay for shipping to and from your location if all instructions regarding the method of packaging and shipping The HIGHERDOSE RED LIGHT NECK ENHANCER are followed.

HIGHERDOSE reserves the right to change the method by which HIGHERDOSE may provide warranty service to you, and your HIGHERDOSE RED LIGHT NECK ENHANCER's eligibility to receive a particular method of service. Service will be

limited to the options available in the country where service is requested. Service options, parts availability, and response times may vary according to country. You may be responsible for shipping and handling charges if The HIGHERDOSE RED LIGHT NECK ENHANCER cannot be serviced in the country it is in. If you seek service in a country that is not the original country of purchase, you will comply with all applicable import and export laws and regulations and be responsible for all customs duties, V.A.T., and other associated taxes and charges. Where international service is available, HIGHERDOSE may repair or replace The HIGHERDOSE RED LIGHT NECK ENHANCER and parts with comparable products and parts that comply with local standards.

LIMITATION OF LIABILITY

EXCEPT AS PROVIDED IN THIS WARRANTY AND TO THE MAXIMUM EXTENT PERMITTED BY LAW, HIGHERDOSE IS NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR CONDITION, OR UNDER ANY OTHER LEGAL THEORY, INCLUDING BUT NOT LIMITED TO LOSS OF USE; LOSS OF REVENUE; LOSS OF ACTUAL OR ANTICIPATED PROFITS (INCLUDING LOSS OF PROFITS ON CONTRACTS); LOSS OF THE USE OF MONEY; LOSS OF ANTICIPATED SAVINGS; LOSS OF BUSINESS; LOSS OF OPPORTUNITY; LOSS OF GOODWILL; LOSS OF REPUTATION; OR ANY INDIRECT OR CONSEQUENTIAL LOSS OR DAMAGE HOWSOEVER CAUSED INCLUDING THE REPLACEMENT OF EQUIPMENT AND PROPERTY. THE FOREGOING LIMITATION SHALL NOT APPLY TO DEATH OR PERSONAL INJURY CLAIMS, OR ANY STATUTORY LIABILITY FOR INTENTIONAL AND GROSS NEGLIGENT ACTS AND/OR OMISSIONS. SOME STATES (COUNTRIES AND PROVINCES) DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

GENERAL

No HIGHERDOSE reseller, agent, or employee is authorized to make any modification, extension, or addition to this Warranty. If any term is held to be illegal or unenforceable, the legality or enforceability of the remaining terms shall not be affected or impaired. This Warranty is governed by and construed under the laws of the country in which The HIGHERDOSE RED LIGHT NECK ENHANCER purchase took place.

HIGHERDOSE or its successor in title is the warrantor under this Warranty.

FCC COMPLIANCE STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference,
- (2) This device must accept any interference received, including interference that may cause undesired operation.

COPYRIGHT © 2024 HIGHERDOSE. ALL RIGHTS RESERVED

HIGHERDOSE
RED LIGHT FACE MASK

HIGHERDOSE.COM

help@higherdose.com

@higherdose