EN、DE、FR、IT、ES

**GLO24K**  
UNLEASH THE POWER OF GOLD

# Skin Rejuvenation Beauty Device User Manual

To Maximize Your Results Go to Back Side

# Skin Rejuvenation Beauty Device

## user manual

### To maximize your results go to back side

GLO24K Skin Rejuvenation Beauty Device is an advanced, triple action beauty device that aims to generate more youthful appearance and assist in rejuvenating and beautifying the skin for a radiant, glowing complexion. This beauty tool is based on cutting-edge Red Light LED, Temperature Controlled Gentle Heating, and Vibration Technologies to help revitalize, revive, and rejuvenate the skin for a flawless, glowing, and youthful appearance. Enjoy Spa-at-Home experience with GLO24K.

#### Triple Action

- LED
- Temperature Controlled Light Heating (40°C)
- Vibration.

#### Remember :

1. Please fully charge the device before use.
2. Please recharge the device if not used for a long time.
3. Full charging time is about 2-3 hours.
4. Auto shut-off is set at 15 minutes for safety reasons.

# Product Features


## LED Skin Rejuvenation

LED (Light Emitting Diode) is a technology that is commonly used in the beauty industry and its benefits are associated with more youthful and radiant appearance.

Stimulates Skin Rejuvenation

Can assist with skin balancing and toning for a more radiant look.
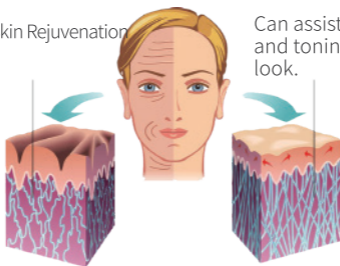
## Temperature Controlled Light Heating & Vibration

GLO24K Skin Rejuvenation Beauty Device, in addition to LED, utilizes Temperature Controlled Light Heating and Vibration / Massage technologies both associated with skin care benefits for a more beautiful, revitalized appearance.

40°C+/-1°C

Controlled Light Heating  
Massage Head

Healthy magnet

Skin Rejuvenation  
Beauty Device

## Your Device

You can charge this device directly on bottom of device or use the base.

## Usage

1. Clean your face, remove make-up, and apply the device only on dry, clean skin.

2. Press A button, Turn on the device, the device will enter the standby mode (blue light).

3. Press the B Button to start up the Heating, Vibration, and LED functions. Red light will show. Once in touch with your skin, these functions are ON and the device head will warm up gradually.

4. Gently massage your face by gliding the device over areas of concern in up/down and circular motions. A recommended use is 10-15 minutes per application 2-3 times a week.

5. Press the B Button to switch to standby mode.

6. To Turn Off the Device press the A button shortly (Not a long press). For best results, apply a moisturizing cream after each application of the device.

Tips: The above methods of care can be adjusted according to your needs.

## Notes & Safety Tips

1. This product is only recommended for use by an adult. Keep it away from children.
2. Wipe off with a dry cloth after each use. Do not wash off or use any cleaning chemicals. This device is not water resistant or waterproof.
3. This device is utilizing a smart sensitive “touch sensor” so it works only when the massage head is in direct touch with the skin and the device is ON.
4. Avoid keeping the device on one skin area for too long and keep on moving it on the skin in circular or upward/downward motions.
5. Close your eyes during application to avoid discomfort.
6. Do not apply this device directly on the eyes. It is designed to be applied on the skin of the face, neck, and around the eyes.
7. Use as needed. Each application for up to 10-15 minutes. 2-3 times a week is recommended.
8. In the rare case of any discomfort or sensitivity stop using this device and contact customer care.
9. If you have any medical condition (such as but not limited to: have a heart pacemaker, epilepsy, pregnancy, suffer from severe medical conditions, etc) please consult with your physician before using this product.

### Product Parameters

Model No: BD1	Material: ABS + Stainless steel
Size: 200*46.5*29.5mm	Net Weight: 240g
Battery Capacity: 1200mAh	Working Voltage: 7.4V
Working Current: <1.5A	Heating Temperature: 40°C+/-1°C
Using Time: 10-15mins/time	Auto shut-off is set at 15 minutes for safety reasons

**Charging:** Flashing indicator lights when charging, blue light when fully charged.

- \* Note: This device is set for an Auto shut-off at 15 minutes for safety reasons.
- \* Charging time may take about 2-3 hours when battery is low.
- \* Flashing lights while ON indicates low battery. Device can not turn ON when the battery is low. Recharge between each use.

**Product Accessories:** Infrared Beauty Device \*1, Charging base \*1  
USB cable\*1, User Manual \*1, cleaning cloth\*1, pouch\*1

## Repair & Maintenance

To clean your device make sure it is turned off and unplugged.

Clean the device with a soft, dry cloth between uses.

This tool is not water proof nor it is water resistant. Do not wash it off with water or any other liquied.

Put the device in a safe place after use. Do not place heavy objects on it.

Store the device in a dry place to avoid moisture.

Recharge device between uses.

If trouble-shooting the device is needed follow the steps below. If the problem is not resolved contact the product service center. Do not disassemble the device.

Breakdown Maintenance	Breakdown Reason	Breakdown Measure
No light therapy or device does not work on the skin	Make sure the device is ON. Make sure the device is touching the skin	Re-adjust gesture to touch skin
Device will flash and can't turn on or does not work	Battery ran out	Please try charging, charging time is about 2-3 hours
Device lights are acting out	Device needs a restart	Turn off the device. Wait a few minutes. Recharge. Turn back ON
Device automatically shuts-off after 15 minutes	Auto shut-off is set at 15 minutes for safety reasons	Wait a few minutes and turn the device back on

**WWW.GLO24K.COM**

Dist. by Y&E Ent., LLC

Atlanta, GA 30338 USA

71-75 Shelton Street, London, UK  
Hauptstrasse 22 Neugersdorf 02727,  
Germany

Warranty:GLO24K warrants all of its beauty tools.  
Check website for details.

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