
Activate your warranty and meet your product

Scan the QR code **within 7 days** of receiving your product to activate your warranty and experience all the product benefits!

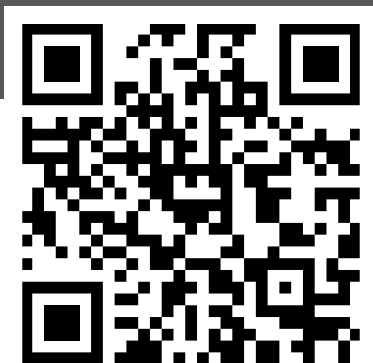Deep Rhythmic Neck Massager With Heat

INSTRUCTION MANUAL AND WARRANTY INFORMATION

NMS-S661H | L-05599, Rev. 1

IMPORTANT SAFETY INSTRUCTIONS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER – TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
- DO NOT reach for an appliance that has fallen into water. Unplug immediately.
- DO NOT use while bathing or in a shower.
- DO NOT place or store appliance where it can fall or be pulled into a tub or sink.
- DO NOT place in or drop into water or other liquid.

WARNING – TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:

- An appliance should NEVER be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual.
- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT sit, lie down or step on the device. This will break the device.
- NEVER drop or insert any object into any opening.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- DO NOT carry this appliance by the supply cord or use cord as a handle.
- To disconnect, turn all controls to the off position, then remove plug from outlet.
- DO NOT use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs, contact Consumer Relations.
- DO NOT use the device in any other way than specified in this manual.
- NEVER operate on a soft surface, such as a bed or couch where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- DO NOT expose a battery pack or appliance to fire or excessive temperatures.
- Exposure to fire or temperature above 80°C may cause explosion.
- Follow all charging instructions and do not charge the battery pack or appliance outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire. Product operating and charging range: 0°C – 40°C.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- DO NOT use outdoors.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by

children or incapacitated persons may be dangerous.

- DO NOT use massager in close proximity to loose clothing or jewelry.
- Keep long hair away from massager while in use.
- DO NOT use a battery pack or appliance that is damaged or modified. Damaged or modified batteries may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.

SAVE THESE INSTRUCTIONS

CAUTION – PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING.

- This product is not intended for medical use.
- Consult your doctor prior to using this product, if:
 - You are pregnant
 - You have a pacemaker
 - You have any concerns regarding your health
- NOT recommended for use by people with diabetes.
- NEVER leave the appliance unattended, especially if children are present.
- NEVER cover the appliance when it is in operation.
- DO NOT use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER use this product directly on swollen or inflamed areas or skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- DO NOT use this product before bed. The massage has a stimulating effect and can delay sleep.
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies.
- This unit should not be used by children or invalids without adult supervision.
- NEVER use this product in automobiles.
- This appliance is intended for household use only.
- This product has an internal, nonreplaceable lithium-ion battery. This battery is not user-serviceable. Please dispose in accordance with local, state, province, and country regulations.
- NEVER use this appliance on or near eyes or other highly sensitive areas.
- All servicing of this unit must be performed by authorized Homedics service personnel only.
- STOP using the device and consult your healthcare provider if you experience an adverse reaction from using the device.
- DO NOT place the device over any open wounds.

INSTRUCTIONS FOR USE

1. Your massager, when opened, may require charging. To charge your massager, plug one side of the USB-C cord into the unit and plug the other end into an adapter. Next, plug the adapter into a 120-volt outlet. Your massager will be fully charged after approximately 4 hours of continuous charging, and the charge will last approximately 3 hours.

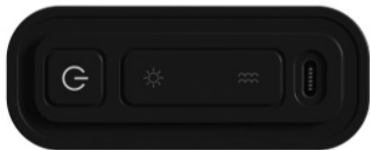
Fig.1

NOTE: This device only comes with a USB-C cord, not an adapter. A USB-C phone charger will also charge the unit. The unit will not operate when it is charging. Use a 5V1A charger, but can not be higher than 5V. Use only with Listed/Certified ITE Power Supply.

2. To turn the massager on, press and hold the power button  for 2 seconds. The LED light indicator will illuminate WHITE.
3. By default, Vibration and Heat is OFF. You can press respective function button once to activate desired functions.
4. Press the vibration button to cycle through 6 vibration modes in different LED indicator. White for Mode 1, Blue for Mode 2, Green for Mode 3, Orange for Mode 4, Purple for Mode 5, Teal for Mode 6. The LED light will be off when vibration is off.
5. To activate soothing heat, press the heat button once. Press again to change toggle low, medium, and high heat. The LED light indicator will illuminate GREEN, BLUE and RED accordingly. To turn heat off, cycle through 3 modes; the 4th press turns off the heat and the LED light will be off.
6. To turn the massager off, press and hold the power button  for 2 seconds. The LED light indicator will turn off.

NOTE: This unit is equipped with a 15-minute auto-off feature for your safety.

LED LIGHT INDICATOR (DURING USE AND CHARGING)

Vibration mode	Press the vibration button to cycle through 6 vibration modes different LED indicator. White for Mode 1, Blue for Mode 2, Green for Mode 3, Orange for Mode 4, Purple for Mode 5, Teal for Mode 6. Turns off when vibration is off.
Heat mode	Heat LED illuminate GREEN (low), BLUE (Medium), RED (High). Turns off when heat is off.
Low power	Power LED flashes WHITE. LED light off when power is depleted.
While charging	Power LED breathes WHITE.
Full charging	Power LED illuminate WHITE. Turns off after unplug from charger.

MAINTENANCE

TO CLEAN


- Be sure to unplug the unit and allow it to cool before cleaning. Use a soft, slightly damp cloth to wipe. Stubborn stains can be removed by wiping gently with 70% alcohol.
- DO NOT use any other cleaners as they may damage the device.
- DO NOT clean the device while it is turned on or charging.
- NEVER allow water or any liquids to come into contact with the product.

TO STORE:

- Store this unit in its box, storage bag, or a cool, dry place.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface.
- Before each use, with the device switched off, clean following the cleaning instructions above.
- DO NOT hang the unit by the power cord.
- 25°C is the most ideal store and operation temperature.

LOCATION OF PARTS AND CONTROLS

Elevate your calm: Massage + heat delivers a full sensory experience

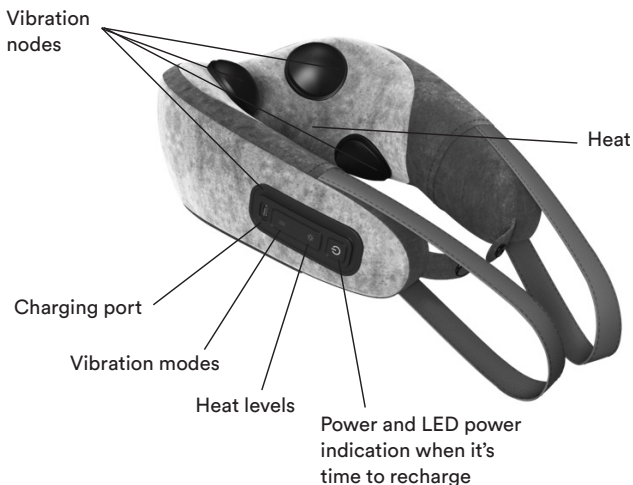SUPPLIER'S DECLARATION OF CONFORMITY

Product Description: DEEP RHYTHMIC NECK MASSAGER WITH HEAT

Model Number: NMS-S661H

Trade Name: Homedics

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

U.S. Contact Information

Company: Homedics, LLC.

Address: 3000 N Pontiac Trail, Commerce Township, MI 48390

8:30am-7:00pm ET Monday-Friday 1-800-466-3342

Homedics is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CAN ICES-003(B)/NMB-003(B)

1-YEAR LIMITED WARRANTY

Homedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of 1 year from the date of original purchase, or the date of delivery, whichever is later, except as noted below. Homedics warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

IMPORTANT: To ensure the most efficient processing of your warranty claim and to stay informed about product updates and new offerings, please register your product at www.Homedics.com/product-registration. Registration is quick, easy, and offers benefits such as faster warranty service, product update notifications, and special promotional offers.

Purchase must have been made from an Authorized Retailer (or their web site) or an Authorized Reseller. Authorized Resellers are Amazon.com, Walmart.com and ebay.com, but only if they say "Shipped and/or Sold by Homedics." **PRODUCTS PURCHASED FROM ANY OTHER ONLINE SELLERS WILL VOID THE WARRANTY.** A copy of a sales or gift receipt is required for all warranty claims. This warranty is valid only for the original purchaser or gift recipient from the date of original purchase, or the date of delivery, whichever is later, and is not transferable.

To obtain warranty service on your Homedics product, contact a Consumer Relations representative by email at cservice@Homedics.com or by phone at 1-800-466-3342 (USA) or 1-888-225-7378 (Canada). Please have the product model number available and be prepared to describe the nature of the problem. Our Consumer Relations representatives are available Monday-Friday 8:30am-7:00pm ET (USA) or 8:30am-5:00pm ET (Canada). Do not send any products to Homedics without contacting us first and obtaining a Return Authorization number. Homedics will, at its sole discretion, repair, replace, or provide a credit for any product that proves to be defective in materials or workmanship during the warranty period. If a replacement product is provided, it will be of equal or greater value and functionality. Replacement may be with a new or refurbished product or components. If the product is no longer available, a product of equal or greater value may be substituted. Any product or part that is replaced becomes the property of Homedics. The warranty period for any replacement product continues from the purchase or delivery date of the original product, not from the date of replacement.

This warranty does not cover: Products purchased from unauthorized retailers, resellers or distributors; damage caused by misuse, abuse, negligence, accidents, or improper maintenance; use for purposes other than those for which the product was designed or intended; damage resulting from improper or unreasonable use or maintenance; failure to follow operating instructions provided with the product; damage resulting from unauthorized repairs, alterations, or modifications; cosmetic damage including, but not limited to, scratches, dents, and broken plastic on ports, unless resulting from defects in materials or workmanship; damage caused by normal wear and tear, including but not limited to, the natural breakdown of colors and materials over extended time and use; damage from external causes such as fire, flooding, earthquakes, or other acts of nature; damage from power surges, power failures, or improper electrical supply; damage to any accessories used in conjunction with a covered product; products with removed, defaced, or otherwise illegible serial numbers or date codes; products sold "AS IS", as "refurbished", or as "reconditioned"; freight, shipping, or transportation costs for returning the product; loss of use or any other incidental, consequential, or indirect damages. Any software included with the product is provided "AS IS" unless expressly provided for in any enclosed software limited warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS OR PROVIDE A WEB COUPON FOR HOMEDICS' WEBSITE IN LIEU OF REPAIR OR REPLACEMENT. TO THE EXTENT PERMITTED BY APPLICABLE LAW, HOMEDICS' LIABILITY SHALL NOT BE GREATER THAN THE PURCHASE VALUE OF THE PRODUCT.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state or country to country. Some states or countries do not allow limitations on implied warranties or the exclusion

or limitation of incidental or consequential damages, so these limitations may not apply to you. For consumers who are covered by consumer protection laws or regulations in their country of purchase or residence, the benefits conferred by this warranty are in addition to all rights and remedies conveyed by such consumer protection laws and regulations. This warranty does not exclude, limit or suspend any rights of consumers arising from non-conformity with a sales contract.

Homedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate Homedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Homedics.

For warranty claims, Homedics will cover the cost of shipping the replacement product for valid claims within the USA. Consumers are responsible for the cost of shipping the product to Homedics.

If you have a dispute, we hope to resolve it as quickly and easily as possible. First discuss your dispute with a Consumer Relations representative. If the Consumer Relations representative is unable to resolve your dispute, you agree that either Homedics or you can initiate arbitration as described in this section (the "Arbitration Agreement"). Arbitration means an impartial third party will hear the dispute between Homedics and you and provide a decision. Binding arbitration means the decision of the arbitrator is final and enforceable. A dispute is any unresolved disagreement between Homedics and you. A dispute also includes a disagreement about this Arbitration Agreement's meaning, application, or enforcement. The arbitrator shall decide any dispute under this Arbitration Agreement. Homedics and you each agree to waive the right to a jury trial or a trial in front of a judge in a public court. This Arbitration Agreement has only one exception regarding the disputes that are arbitrable: Either Homedics or you may still take any dispute to small claims court. This Arbitration Agreement is governed by the Federal Arbitration Act. The American Arbitration Association (AAA) will administer each arbitration and the selection of arbitrations will be according to applicable AAA rules. The arbitration will be held in the state of your residence. Neither Homedics nor you will be entitled to join or consolidate disputes by or against others as a representative or member of a class, to obtain relief in any arbitration in the interests of the general public, or to act as a private attorney general (the "Class Action Waiver"). Homedics will cover the cost of arbitration for any non-frivolous claims, but not the claimant's attorney fees. If any provision related to this arbitration agreement is found to be illegal or unenforceable, then such provision shall be severed from the arbitration agreement, but the rest of the agreement shall remain enforceable and in full effect. Homedics or you each can exercise any lawful rights or use other available remedies to: Preserve or obtain possession of property; exercise self-help remedies including setoff rights; or obtain injunctive relief (including public injunctive relief), attachment, garnishment, or appointment of a receiver by a court of competent jurisdiction. The substance of any disputes where public injunctive relief is available shall be decided by the arbitrator. Only if the claimant succeeds on its claim permitting the remedy of a public injunction may such claimant request that a court of competent jurisdiction enter an injunction in conformity with the arbitral award. You have the right to opt out of this Arbitration Agreement. If you opt out, neither you nor Homedics can require the other to participate in an arbitration proceeding. To opt out, you must contact Homedics by email at the address specified below no later than 30 calendar days from the date of the first consumer purchaser's purchase of the product. You must provide: (a) your name and address; (b) the date on which the product was purchased; (c) the product model number; and (d) the serial number found on the product. Opting out of the Arbitration Agreement shall not be construed as a modification of any other terms, including the Class Action Waiver.

For more information regarding our product line in the USA, please visit: www.Homedics.com. For Canada, please visit: www.Homedics.ca.

FOR SERVICE IN THE USA

cservice@Homedics.com

8:30am-7:00pm ET Monday-Friday

1-800-466-3342

FOR SERVICE IN CANADA

cservice@Homedics.ca

8:30am-5:00pm ET Monday-Friday

1-888-225-7378

©2025 Homedics, LLC. All rights reserved.

Homedics is a registered trademark of Homedics, LLC.

Distributed by Homedics, LLC, 3000 N Pontiac Trail, Commerce Township, MI 48390

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©2025 Homedics, LLC. Todos los derechos reservados.
Homedics es una marca registrada de Homedics, LLC.
Distribuido por Homedics, LLC, 3000 N Pontiac Trail, Commerce Township, MI 48390