

CONGRATULATIONS ON YOUR NEW AIR COMPRESSION BOOTS FLEX.

Please read all instructions and warnings prior to use. Keep this manual as it contains important information for future reference.

SHARPER IMAGE® NAME AND LOGOS (INCLUDING THE 77 LOGO) ARE TRADEMARKS OR REGISTERED TRADEMARKS. DISTRIBUTED BY AND © 2024 MERCHSOURCE, LLC. IRVINE, CA 92618. ALL RIGHTS RESERVED.

For customer support, please visit <https://msccustomerca.re.thethreesixtygroup.com/> or call 1-800-374-2744.

MODEL NO.: 1019355, 1019356, 1019357

Printed in Vietnam.



SHARPER IMAGE®

AIR COMPRESSION BOOTS FLEX

ITM. / ART. 2816330

CARE INSTRUCTIONS

- Spot clean only. Unplug and wipe with a damp cloth.
- Air-dry only. Never place in a dryer.
- Allow to dry fully before connecting to power sources.
- Do not submerge in water.
- Do not use detergent or soap when cleaning.
- Hang the boots or lay them flat for long-term storage.
- Never crumple or wad the boots.
- Do not crimp or forcefully bend any of the included hoses.

SAFETY WARNINGS

- Not a children's item. Keep out of the reach of children.
- As with any massage, bruising may occur. Test your comfort level and check targeted areas often. Stop your massage immediately if you feel any increase in physical discomfort or pain.
- If you experience numbness, tingling, or loss of sensation, stop your massage immediately. Seek professional medical advice if these symptoms persist.
- Do not fall asleep while using the massager.
- Do not use attachments of any kind other than those that may be included.
- Do not leave this massager unattended when it is plugged into a power source.
- Do not expose to extreme temperatures or moisture.
- Keep away from sources of water and heated surfaces.
- Never touch the power cable or ports with wet hands.
- Do not carry this item by the power cable, hoses, or attachments.

- For indoor use only.
- Thoroughly inspect product before every use. If the product appears worn, frayed, splintered, cracked, or broken in any manner, discontinue use and discard immediately.
- This product is not intended to treat, cure, or prevent medical conditions.
- Do not use this item if you are pregnant or if have an injury or illness. Consult your physician if you have a medical condition that may prevent safe use of this item.

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Terms of Warranty

- 2 Year Limited Warranty : Device and Charger.
- 90 Days: All other accessories.
- One-Year Limited Warranty: The original unit and/or non-wearable parts deemed defective. It's at the sole discretion of MerchSource to determine warranty application.
- In the event replacement unit is issued – the warranty coverage does not extend, the original date of purchase will remain as starting point for the one-year policy.
- Use of non-original battery or charging cable is not allowed. Any use of any non-DEM components will void warranty and any responsibilities.
- Product must be purchased at approved retailer. No warranty for renewed/refurbished units.

Guidelines for all Warranty-Returns:


1. Proof of Purchase required from approved retailers.
2. E-Boy, Coffer-up, Discount Stores (Dollar Store type) not approved, and no warranty/returns can be applied to our products if purchased from such retailers-resellers.
3. Customer shall provide image of original receipt matching sku- product details.
4. No financial reimbursement – only product replacement. If originally purchased product is out of stock or no longer available, MerchSource reserves the right to make product substitution to remedy. No refund is allowed.
5. All product shipments are shipped using Ground Service – lowest cost available.
6. Any returned products must be received within 14 days.
7. It's always recommended to consumer to return the product to the original place of purchase.
8. Consumer information is required and all product details including batch numbers. Consumer information is not shared or used for any marketing.

Returns Process-Guidelines:

- If product has failed due to potential safety, potential injury – such product must be returned to MerchSource.
- In most cases defective products/components/accessories do not need to be returned but consumer must wash responsibly and/or mutilate product to prevent further use.
- Product Mutation (1): Consumer shall provide current image of product before and after mutation. If original cable is supplied – consumer shall disconnect charger and cut the cord. Image is required showing proof.
- Product Mutation (2): Consumer is asked to trash product and all its accessories (image is required of product in trash displaying proof of such action).

WHAT YOU GET

- a. Air Compression Boots Flex (x2 boots)
- b. Control Unit
- c. Air Hose Extension
- d. Power Adapter with DC Jack

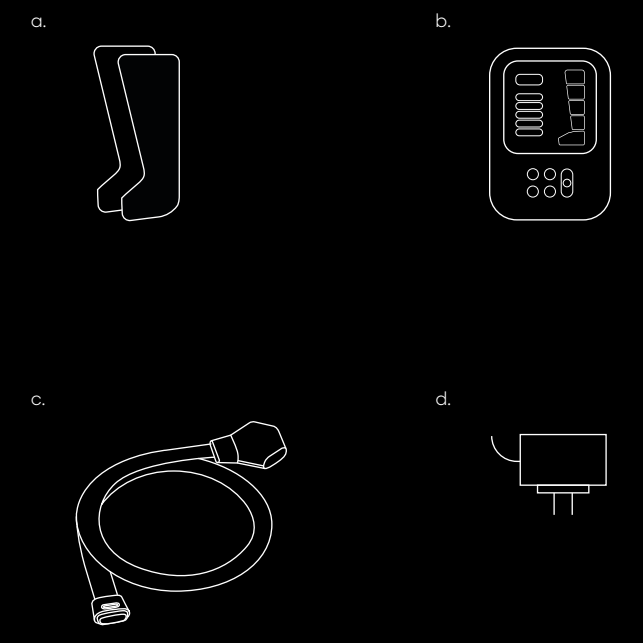


THE FUTURE OF RECOVERY

01 GETTING STARTED

1. Plug the barrel connector on the power adapter cord into the power input port on the control unit.
2. Plug the air hose extension into the control unit.
3. Plug the power cord to a wall outlet.

Note: Do not use the power adapter to charge or power other devices.

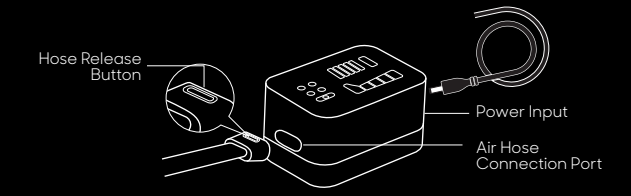

4. Unzip the boot or boots you want to use for your massage.
5. The boots are labeled L (left leg) and R (right leg).
6. Find a comfortable place to sit or recline with your legs extended for the duration of your massage and remove any footwear.
 - Wearing trousers and socks are recommended.
7. Pull each boot over left & right legs and zip the boot up all the way to the top of the zipper track.

IMPORTANT:

- When using only 3 or 4 zones manually, do not roll over the boot or zip it mid-way to prevent damage.
 - Never unzip the boots when they are inflated, as this will damage them.
 - Only zip or unzip the boots when they are fully deflated.
 - Attempting to operate the zipper when not deflated can cause damage.
- Lock the zippers in place by flipping the zipper tabs down fully flat against the zipper track. The tabs should point towards the toes of the boots. Failure to properly lock the zippers may cause the boot to open during inflation, reducing compression.

IMPORTANT: Never attempt to stand or walk while wearing the boots. The boots are intended for massage only. They are not designed for stability or mobility. Using the boots for purposes other than their intended.

02 CONNECTING THE AIR HOSES

IMPORTANT: Connecting the Air Hoses

1. Orientation:

- Ensure the hose release button is facing upwards when connecting each hose.
- Press each connector into the connection port slowly and firmly until it clicks securely into place.

2. Using Both Boots:


- Plug the single connector end of the air hose extension into the control unit's air hose connection port.
- Connect both boot air hoses to the connection ports on the other end of the air hose extension.

3. Using a Single Boot:

- Option A:
- Connect the boot air hose directly to the control unit.


Option B:

- Use the air hose extension:
- Connect the boot air hose to one port of the extension.
- Seal the unused connection port with the rubber plug on the bottom of the connector junction (the top side has the Sharper Image logo).
- Connect the single connector end of the extension to the control unit.

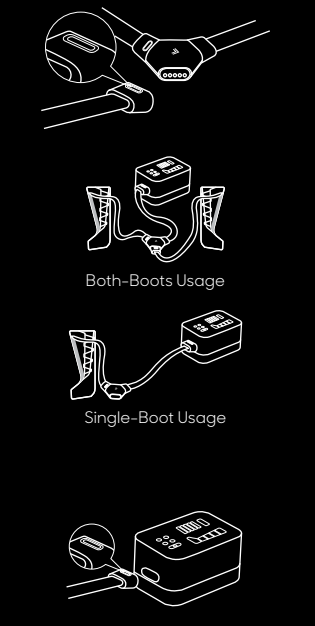

03 DISCONNECTING AIR HOSES

- Press the hose release button while gently pulling the connector away from the connection port.
- Always grip the plastic connectors directly.
- Never pull or push using the hoses themselves.

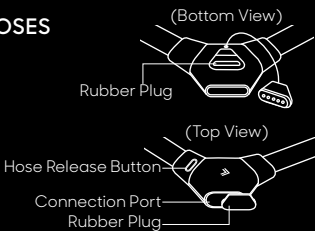

04 OPERATING THE MASSAGER:

1. Power On/Off:

To turn on: Press and hold the power button.

To turn off: Press and hold the power button again.

2. Timer Settings:

After compression starts, short press the timer button to adjust massage duration.

Each short press cycles through +10 minutes, +20 minutes, return to initial time.

3. Massage Modes:

Press the mode button to cycle through 5 massage modes:

- a. Pulse (default): All zones compress and decompress.
- b. Gradual: Zones compress sequentially, then all decompress.
- c. Sequence: Zones compress and decompress sequentially.
- d. Zones: Select up to 5 zones for targeted pulse massage.
- e. Manual: Select up to 5 zones and manually adjust compression.

4. Zone Selection (for Zones and Manual modes):

- a. Use up, down, and select buttons to choose compression areas.
- b. Press and hold the select button to start the massage.

5. Compression Levels:

Quick press the compression button to cycle through 5 compression levels.

6. Ending the Massage:

- a. Quick press the power button during massage & allow the boots to decompress gradually.
- b. Detach the hose from the control box and disconnect the boots from the extension connector
- c. Gently clean the boots with isopropyl alcohol / alcohol wipes
- d. Hang the boots to air dry.

