LED Light Therapy Device

User Manual

Please read this manual carefully before using the product

COMPANY INFORMATION

At Numiere, we believe skincare is a journey, not about fixing flaws. Lines, spots, and texture are part of your story, and we celebrate them. With years of expertise in the skincare and beauty industries, we're redefining who has access to the best tools and technologies to support their unique journey.

Our mission is to make high-quality skincare accessible to everyone. Numiere's tools are simple, effective, and designed for use anywhere, by anyone. We're here to empower you with restorative rituals and advanced innovations that enhance your skin's natural renewal process—wherever you are in your journey.

GENERAL DEVICE INFORMATION

The Numiere Feather LED Decollete Mask is an medical device that uses light therapy to effectively treat wrinkles around the decollete.

DESCRIPTION OF THE DEVICE

The Numiere Feather LED Decollete Mask combines two powerful light therapy treatments in one innovative device:

1. Red Light Mode – Uses red LED light, often selected for its potential effects on the appearance of fine lines and skin tone.
2. Near-Infrared Light Mode – Uses near-infrared wavelengths designed to reach deeper layers of the skin.

Energy:	Max 70mw/cm2
LEDs:	58pcs
LED qunntity:	Deep red:630nm+/-10nm
	NIR (Near-Infrared):830nm +/-10nm
Program Time:	10/15/20 minutes
Power:	Max3.5w
Battery Capacity:	1200mAh
Operating Conditions:	Operating Temperature: 0°C to 40°C
Storage and TransportTemperature:	-5°C to 55°without condensation
Coperating Humidity:	30% to 95%
Storage and Transport Humidity:	20% to 85% without condensation
Charging supply:	5V ==1A
Charging cable :	TYPE C
Charging cable length :	80cm
Charging time :	about 90 minutes
Time of use :	80minutes
Recommended use :	1 time per day,3-7times per week
Materials :	Silicone ,ABS+PC,ABS

DEVICE ACCESSORIES

1. Feather LED Decollete Mask&Controller
2. Adjustable, Detachable neck Strap
3. Manual
4. Charging Cable

x1

x1


x1

x1

Important Features:

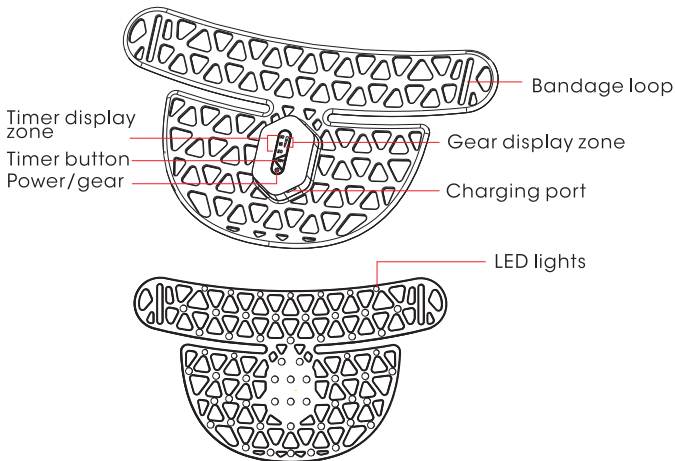ADVERSE REACTION

If you have any adverse reaction from using this device, such as discomfort, skin discoloration, skin irritation, or other concerns, stop using it immediately and consult your health professional. Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your health professional.

If you experience a persistent headache and/or other persistent painful sensations while using this device, stop using it and consult with your health professional

CONTRAINDICATIONS

- DO NOT use the device if you have any active cancer or cancerous lesion.
- DO NOT use the device if you are pregnant, planning to become pregnant, or breastfeeding. The device has NOT been tested on pregnant or breast-feeding women and therefore the risk to the fetus, neonate or pregnant women is unknown.
- DO NOT use the device if you are under the age of 22.
- DO NOT use the device over areas with phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to silicone or plastics. The device contains plastics and silicones that come into contact with the skin during use. When in doubt of allergies, contact your health professional before using the device.
- DO NOT use the device if you suffer from any photosensitive disorder (sensitization to light)
- DO NOT use the device if you suffer from any serious skin condition such as severe acne.

- DO NOT use the device if you are taking any medication or drugs that can cause photosensitivity. If you are not sure about your medication, contact your health professional.
- DO NOT use the device if you suffer from light induced headaches.
- DO Not use the device if you have a medical history of seizures or epilepsy triggered by light.
- DO NOT use the device if you have the Fitzpatrick Skin in Type V and Type VI.
- DO NOT continuously press the ON/OFF button.

Note: We highly suggest being cautious when taking any (skin-related) medications. If in doubt, consult your health professional.

Note: If you are under the care of a physician, consult them before using this device.

WARNINGS

- STOP using the device and consult your health professional if you experience an adverse reaction from using the device.
- STOP using the device if you experience any discomfort or difficulty breathing or swallowing occurring from application of the device on the neck.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness) while using the device.
- STOP using the device if its overheating. Turn it off immediately.
- STOP your treatment if it feels hot or your skin feels uncomfortably warm.
- DO NOT let children play with the device. Cleaning and

charging the device shall not be done by children.

- DO NOT use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- DO NOT use the device on spots without pigmentation or areas on the skin that have vitiligo.
- DO NOT shine the light from the LEDs directly into your eyes or someone else's eyes. This may cause temporary loss of sight.
- DO NOT place the device on any open wounds.
- DO NOT use the device outdoors to prevent damage to the product.
- DO NOT sit, lay down on, or step on the device. This will break the device.
- DO NOT drop or insert any object into any openings on the device.
- DO NOT use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs, contact customer support.
- DO NOT use or store the device near hot surfaces.
- DO NOT charge the device when it is close to or might come into contact with water.
- DO NOT submerge the device in water.
- DO NOT store the device in extreme heat or cold.
- DO NOT place the device in direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. It contains batteries that are non-replaceable.
- DO NOT perform maintenance on or take apart the product. In case of malfunctioning or damage, contact customer support.
- DO NOT use the device if it is damaged or if the charging

cable is damaged.

- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than specified in this manual.
- DO NOT use the device if you have had a recent surgery or other surgical procedures.
- DO NOT use the device over a bacterial infection, as this may result in the worsening of symptoms.
- DO NOT use the device on painful areas. Consult your health professional when you encounter painful areas during use.
- DO NOT use the device on or near children. Do not let it be used by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned, or stored near children or those with reduced physical or mental abilities.
- DO NOT use the device with any USB cable other than the one provided by Numiere.
- DO NOT use the device over areas of skin that lack normal sensation. Only use the device on normal, intact, clean, and healthy skin.
- DO NOT use the device while driving, operating machinery, or during an activity in which voluntary muscle contractions may put users or the public at any risk of injury.
- DO NOT use the device on swollen, infected, or inflamed areas or skin eruptions.
- DO NOT use the device in high humidity locations such as the shower or sauna.
- DO NOT exceed the recommended treatment time when using the device.

- If there has been sunburn on your skin in the past two weeks, it is recommended not to use this product.
- This device is intended to be used by one person, we recommend to avoid sharing with others.

Note: The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cable is damaged, discontinue use. Contact Numiere customer support to arrange a cable replacement.

PRECAUTIONS

- Do not scratch or damage the light-emitting surface.
- Not operating the device in accordance with the instructions in the manual will void your warranty.
- Read all instructions carefully before using this device.
- Store the device according to the protocol outlined in this manual.
- This device is intended to be used by one person only.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your health professional if you are unsure about your health status.
- The device's battery is not removable. If you are disposing of the device, do not burn it or throw it in the trash. Dispose of the device according to local regulations.
- Inspect the device and do not use it if you see any signs of damage, such as internal battery leaks. If anything seems out of the ordinary, contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cord is damaged, discontinue use. Contact customer support to arrange a replacement.
- Consult your health professional before using the device

if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anesthetic, or desensitized areas.

- We highly suggest using caution if you are taking any medications. If in doubt, consult your health professional. If using the device causes persistent pain or discomfort, discontinue use immediately.
- There may be a perception of flashing lights during the use of the device, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your device, consult your health professional.

DISCLAIMER

This device is not intended to diagnose, cure, or prevent disease. Numiere makes no claims, representations, or warranties regarding the ability of this device to cure any physical, skin, or mental condition. A qualified health professional should be consulted regarding any condition requiring medical attention.

RISKS AND BENEFITS

POTENTIAL RISKS

- Immediately after treatment, skin may appear red. Skin should return to normal within 1-2 hours.
- Temporary headaches may develop immediately following the light treatment. Headaches can be treated with an appropriate pain killer and rehydration.

POTENTIAL BENEFITS

- Red Light Therapy: reduces the appearance of fine lines and wrinkles.
- Near-Infrared Light Therapy: helps manage suboptimal

local blood flow and circulation while addressing muscle and joint pain

LIGHT SENSITIVITY TEST

We highly recommend users perform a light sensitivity test before their first use of the device.

1. Place the device on the inner forearm.
2. Turn on the device to start a treatment.
3. Wait until the treatment is done. The full treatment will take 3 minutes to complete.
4. If your skin feels hot or uncomfortable during the treatment, stop immediately! You might be sensitive to light.
5. The device will automatically turn off after the treatment is done. Remove the device and wait for at least 24 hours before your next treatment.

Check your skin regularly over those 24 hours. If your skin appears red and/or itchy, the device is unsuitable for you. In case of sensitivity to the light, please contact customer support.

HOW TO USE THE DEVICE

GENERAL USE

- a) Hold the ON/OFF button 2s to turn on the device. The indicator ring around the controller will turn on white. short press the ON/OFF button to adjust the light level;
- b) Press the timer button to choose different timer.

TREATMENT PROTOCOL

Use the device for 3 minutes per day, 3 days per week for a minimum of 4 weeks. After that, you can continue using

the device as desired.

Always be mindful of how your skin reacts to beauty devices, and decrease usage if redness or irritation occur.

CHARGING

Charge your device fully up to 2 hours before your first use. Only use a 5V adapter or USB port (avoid fast chargers as they can damage the battery).

To charge the device, attach the charger to the device. Insert the USB-Aside of the cable into a USB power supply. The LED indicator on the device will start to blink to indicate charging. When the device is fully charged, the LED indicator will stay on.

If the LED indicator remains on when you first connect the charging cable, then your device is already fully charged and you may begin using. If the light is blinking please charge until the light remains steady.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use. The device will not work while being charged.

TROUBLESHOOTING

If the following suggestions do not fix the problem or the problem reoccurs, and your device is within the warranty period, please return it to Numiere.

- DO NOT attempt to investigate any further.
- DO NOT attempt to modify the device in any way.

The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged,charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged,because the cable was not properly connected to the device.	Connect the cable properly to the device and charge the device.
The device doesn't switch on.	The power outlet didn't work during charging.	Try a different power outlet.
	The ON/OFF button has not been pushed for 2 full seconds.	Push the ON/OFF button for 2 full seconds.
The device switches off unexpectedly.	The battery level is too low.	Charge the device.

CARE AND STORAGE

The lifetime of the lithium battery can be maximized by allowing it to fully charge and discharge during regular use. Allow the device to deplete before charging again. Do not charge after every session.

- Store the device at room temperature in dry conditions.
- Store the device in the supplied bag. Make sure the device will not bend, stretch, or be crushed while being stored or when traveling.
- Avoid storing the device in dusty conditions, direct sunlight, or damp and humid conditions.
- This device contains no serviceable parts.
- Before each use, ensure the device is switched off and clean it, following the cleaning instructions below.

IMPORTANT

- The product has a limited lifecycle of two years under nominal use. Following the treatment protocol for two years is considered nominal use.
- The device must be handled with care.

CLEANING

- Clean the device with a water-based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with a 70% alcohol wipe.
- DO NOT use any other cleansing devices as these may damage the device.
- DO NOT clean the device while it is turned on.
- DO NOT clean the device while it is charging.

DISPOSAL

- DO NOT dispose of the device with household waste. The device should be recycled. The device contains a lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

WARRANTY

Limited Warranty

Numiere warrants that this product is free of defects in materials and workmanship for one (1) year from the date of purchase. This warranty is valid only in accordance with the conditions set forth below:

- Normal wear and tear is not covered by this warranty.
- This one (1) year warranty applies to consumer use only and is void when the product is used in a commercial or institutional setting.
- This warranty extends only to the original purchaser and is not transferable. In addition, proof of purchase must be demonstrated. This warranty is void if the product has been subject to accident, misuse, abuse, improper maintenance or repair, unauthorized modification, or any other use not found in the printed instructions.

This limited warranty is the only written or expressed warranty given by Numiere. Any implied warranties on the product (including but not limited to any implied warranties of merchantability or fitness for a particular purpose) are limited in duration to the duration of this warranty. Some statutes do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you. Repair or replacement of the product is the exclusive remedy of the consumer under this warranty. Numiere shall not be liable for any incidental

or consequential damages for breach of this warranty or any implied warranty on this product. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Procedure for Warranty Repairs or Replacement:

If warranty service is necessary, the original purchaser must contact Numiere for warranty service instructions at support@Numiere.com

For your protection, it is suggested you insure the shipment of the device back to Numiere. Damage that occurs during shipment is not covered. Packages that are missing proper documentation and accessories may not be processed.

SAFETY

The device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The device has been tested to or complies with the following standards:

- EMC and Low voltage Directive IEC 60601-1, -1-2, -1-6, -1-11, -2-57.
- RoHS Directive 2011/65/EU.
- IEC60601-1-6 and IEC62366 usability of Medical Device.
- IEC 64271 Photobiological safety of lamps and lamp systems.

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