



HEAT SEALED POWDER COAT STEEL 14 SCOPE DRYING CABINET, E-LOCK

SC8036DREDP-14

Scope Cabinet Features:

- Meets or exceeds Joint commission, AORN, SGNA, and AAMI guidelines.
- Positive pressure fully redundant fan system
- Industries only unit with 3 intake and 2 exhaust fans
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Industry first replaceable HEPA Filter system without removal of scopes
- Holds 16, 14, 12 or 8 scopes with secure adjustable mounting system
- Medical grade electrical system
- Tempered glass double doors
- Easy to clean removable drip tray

The Sure-Dry Scope Cabinet includes the industry best Positive Pressure HEPA Filtered Drying system. Inside air is replaced inside the cabinet 4 times per minute.

Scopes
not included

*Optional dimmable LED
lighting also available!*

Simple “Plug and Play” upgradeable with Dri-Scope Aid®, industry leading internal lumen drying system


HARLOFF
CLINICAL STORAGE BY CAPSA HEALTHCARE

Engineered &
Manufactured
in the U.S.



View similar products online at
<https://www.harloff.com/suredry-scope-cabinets/>

Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300

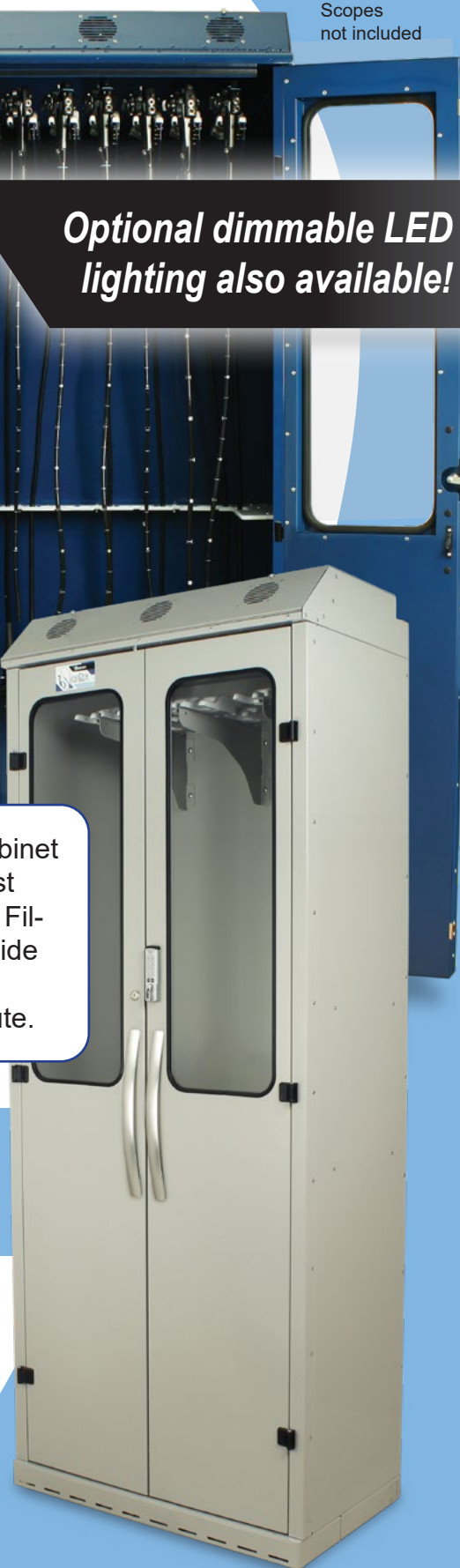

STAINLESS STEEL 14 SCOPE CAPACITY CABINET, E-LOCK



SureDry
SCOPE CABINETS

SCSS8036DREDP-14

Scopes
not included

Scope Cabinet Features (continued):

- Multiple locking options available including badge entry
- Scope cabinets are delivered fully assembled
- Harloff has Industry Best 12 year manufacturers warranty for Sales to US and Canada
- Overall cabinet dimensions:
93"H x 36"W x 24"D

Color Options:

- Stainless Steel
- Heat sealed powder coat low carbon steel, available in over twenty color options
- View colors online:
<https://www.harloff.com/finish-options/>

HARLOFF
CLINICAL STORAGE BY CAPSA HEALTHCARE

Engineered &
Manufactured
in the U.S.



View similar products online at
<https://www.harloff.com/suredry-scope-cabinets/>

**Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300**



Timer:

- Timer sets drying time with both doors closed, and resets when doors open.




- Default is 90 minute runtime, but can be set from 10 minutes to 180 hours.

Why is a positive pressure system important?

The unique design of the SureDry Scope Cabinet includes both top fans for intake and lower fans for exhaust. This push/pull design creates a vortex of air that flows through the cabinet ensuring maximum drying efficiency.

Intake: Intake airflow comes through separate HEPA filters (which remove 99.7% of particulates, .03 Microns) driven by fully redundant fans producing a positive pressure vortex from the cabinet top. This creates an even circular motion of airflow over all scopes.

Exhaust: The lower exhaust fans pull the circulated air out of the cabinet ensuring the flow of air is consistent from top to bottom, providing a complete drying system for the entire scope. This push/pull multiple fan system (an industry first) ensures a positive pressure system to meet all major guidelines and it ensures forced drying and equal drying levels across the full length of each and every scope.

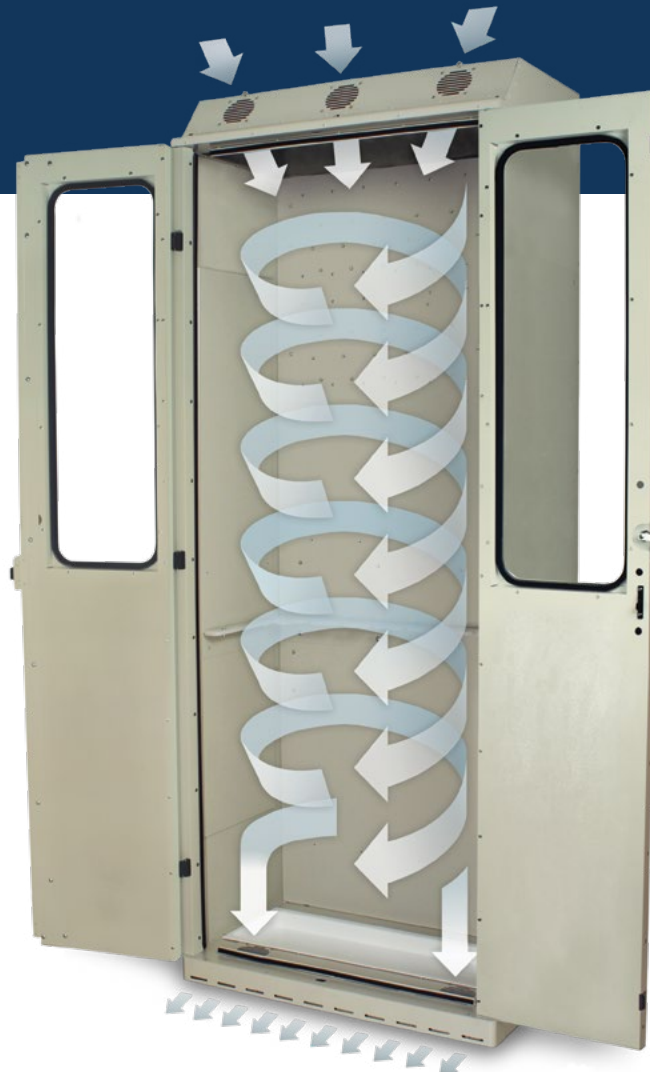

ST91:2021

SureDry Cabinets are a Smart Choice to Meet and Exceed ANSI/AAMI ST91:2021 and Other Guidelines

The references below are taken from these guidelines. Channel drying available only with the addition of the Dri-Scope Aid® system.

☒ **Scopes Stored Vertically**

11.2.1 "stored vertically, the insertion tube should be as straight as possible, with the distal tip hanging freely."
11.2.1 "When flexible and semi-rigid endoscopes are hung in the vertical position, coiling or kinking is prevented."

☒ **Hanging Freely Without Coiling**

11.2.2 "should be stored in a cabinet that is of sufficient height, width, and depth to allow flexible endoscopes to hang vertically without coiling and without any components touching the bottom of the cabinet."

☒ **Channel Drying + Forced Air**

8.2.5.1 "Flexible endoscopes with channels should be dried for a minimum of 10-minutes with pressure-regulated forced instrument air or a minimum of HEPA-filtered air." | 8.2.5.2 Use pressure-regulated instrument air or HEPA-filtered air to dry the channels in accordance with the manufacturer's written IFU.

☒ **HEPA Filtered Positive Pressure**

11.2.2.2 Cabinets that enable circulation of HEPA-filtered or instrument air through the cabinet at... positive pressure.
NOTE: SureDry cabinets feature an industry leading 3 intake and 2 exhaust fan system.

☒ **Easy Scope Identification***

11.2.3 "Protocols should be developed to help ensure that users can readily identify an endoscope that has been processed and is ready for patient use." * Available with Tracker-equipped SureDry scope cabinets

☒ **Single-Use Connectors**

Use of single use connectors ensures no cross-contamination as well as simplifies staff efficiency.

☒ **Clinical Documentation**

ASGE "Drying the interior of an endoscope should use a sufficiently prolonged flow of medical air through all accessible channels, ideally simultaneously and for at least 10 minutes, for greatest efficiency." Per the study, "Multisociety guideline on reprocessing flexible GI endoscopes: 2016 update."
NOTE: SureDry cabinets are the only clinically proven vertical drying cabinet.

10 Reasons to Choose **HARLOFF** for Your Scope Cabinets

CLINICAL STORAGE BY CAPSA HEALTHCARE



Industry Best Positive Pressure HEPA Filtered Drying System. Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

1

A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines

Proudly Made in the USA



Harloff products are proudly **Designed, Engineered and Manufactured in the United States** at our Colorado Springs, CO facility.

6


Dri-Scope Aid® Fully Retrofittable in Field

All SureDry scope cabinets are fully designed to be field installed with the industry best Dri-Scope Aid® lumen drying system, which provides HEPA filtered air to each internal scope channel.


2

Easy Access to HEPA Filters

Change HEPA filters **without opening the cabinet.** Save time and ensure no cycle interruptions.

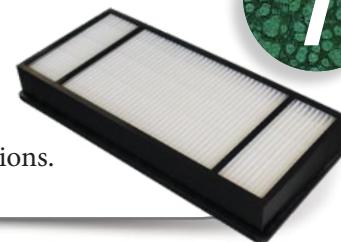

7

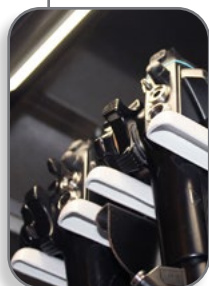

LED Lighting Option

Brighten your cabinet and easily inspect scopes with fully dimmable optional overhead LED lighting which triggers on door open/close.

3

Keyed or Keyless Locks

Options include key locks, basic electronic push button locks, advanced keypads and proximity readers. View locking options online at www.harloff.com/locking-options/



8

Hooks and Holders for Specialized Scopes

Major brands to include **Fuji, Olympus, Pentax, Karl Storz and others.**



4

Scope Cabinet Finish Options

Choose Stainless Steel, Heat Sealed Powder Coat in High Gloss or Textured Low Gloss, or upgrade to an **antimicrobial finish.**

9

12 Year Warranty

Our Manufacturer's Warranty is Unmatched in the industry. We build our cabinets to last and we stand behind them.

5



Superior Color Choices

Over 20 Colors to choose from. Cabinets have a heat sealed powder coat finish inside and out. View color options online at www.harloff.com/finish-options/

10