



TRIO.BASTM

BIOLOGICAL AIR SAMPLER

**The most complete line
of active microbial air
sampling instruments
available**



TABLE OF CONTENTS


| | | | |
|---------|--------------------------------|---------|---------------------------------|
| 1..... | Origins | 20..... | TRI CLAMP Satellite |
| 2..... | Principles | 21..... | TRIO.GAS |
| 3..... | Regulations | 22..... | FLO.GAS and VERI.GAS |
| 5..... | MINI | 23..... | Calibration Control Equipment |
| 6..... | MINI for CompactDry™ | 24..... | DAILY SHIFT HEADS |
| 7..... | MONO | 25..... | Data Transfer Options |
| 8..... | DUO | 26..... | Biological Air Sampler Software |
| | MONO/DUO | 27..... | Data Transfer Software |
| 9..... | Data Transfer Cable Models | 28..... | Comparing AS & BAS Software |
| 10..... | AIRBIO ONE and DUO | 29..... | CFU Photo Camera |
| 12..... | AIRBIO ONE RAPID-VIRUS | 32..... | Accessories |
| 14..... | Cleanrooms Monitoring Question | 32..... | Multi Holder Carts |
| | RABS ISOLATOR and | 34..... | TRIO.SETTLE |
| 16..... | MULTIFLEX Systems | | Prepared Culture Media |
| 18..... | Remote Aspirating Tube System | 35..... | and CompactDry™ |
| 19..... | Remote Aspirating Head System | | |

Applicable Industries:

- Pharmaceutical/Biotechnology
- Compounding Pharmacy
- Medical Device Manufacturer
- Sterile Product Packager
- Hospital/Surgical
- Food/Beverage/Dairy
- Cosmetics/Personal Care
- Cannabis
- Tissue Lab
- Aerospace




The Discovery of Airborne Organisms

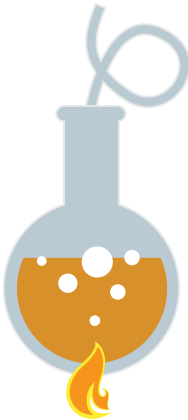


BROTH IS BOILED



**BROTH REMAINS FREE OF
MICROORGANISMS**

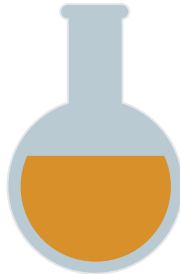

**CURVED NECK
IS REMOVED**

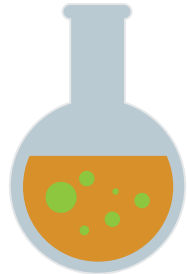

**MICROORGANISMS
GROW IN BROTH**

Lazzaro Spallanzani in the 1700's, and Louis Pasteur in the 1800's, were the two scientists who first demonstrated the presence of microorganisms in the air after several years of experimentation. After three centuries, it is now possible to perform the same test in a few minutes. Presenting the new generation of microbial air samplers, by the original engineers of the first portable air samplers introduced over 45 years ago.

Spallanzani



Pasteur



Principles


Principle of microbial air sampling by impact method on agar plate

Air, containing microbe-carrying particles, is aspirated and accelerated through an intake and directly towards the surface of a plate. As the air turns away from the agar surface, the microbe-carrying particles that cannot follow the flow are impacted. The plate containing agar is then incubated for the appropriate time and at the proper temperature. The resulting Colony Forming Units (CFU) are counted to evaluate the number of microbe-containing particles collected from a specific volume of air.

How the microbe-carrying particles impact on an agar surface:

The aspirated air passes through an intake of the sampler head at a velocity of " U " and, as it approaches the agar surface, it turns. The arc of the turning circle has a radius of " r " which is assumed to be the same as the radius of the intake nozzle. The velocity around the curve is assumed to be " U ".

The microbe-carrying particle travels along the streamline and experiences a centrifugal force that causes it to move toward the agar surface of the plate.

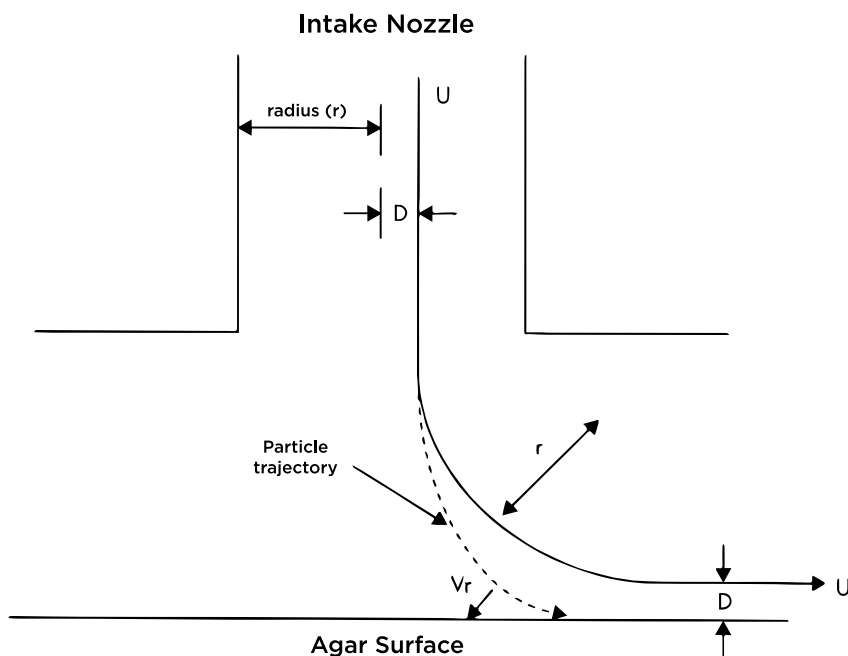

(Fig. 1) Impaction of a particle on a surface after exiting a nozzle

Regulations

FDA Guidance:

According to FDA Guidance for Industry *Sterile Drug Products Produced by Aseptic Processing- Current Good Manufacturing Practice* on Active Microbial Air Monitoring: "...Assessing microbial quality of air should involve the use of active devices... We recommend that such devices be used during each production shift to evaluate aseptic processing areas at carefully chosen locations. Manufacturers should be aware of a device's air monitoring capabilities, and the air sampler should be evaluated for its suitability for use in an aseptic environment based on collection efficiency, cleanability, ability to be sterilized, and disruption of unidirectional airflow."

USP Regulation <797>:

Issued by the non-profit US Pharmacopoeia (USP) and endorsed by the Joint Commission on Accreditation of Healthcare Organizations (JCAHO), USP Regulation <797> is the first enforceable standard for sterile compounding. Originally enacted on January 1, 2004, with the latest version effective November 1st, 2023.

USP <797> is a broad regulation covering a variety of pharmacy policies and procedures designed to reduce the number of patient infections due to contaminated pharmaceutical preparations. USP <797> contains specific requirements for ongoing air and surface evaluation to ensure product sterility and safety for compounded sterile preparations (CSPs).

Recommended Action Levels for microbial contamination:

Viable Air Sampling

| Classification | CFU/1000 liters of air |
|----------------------|------------------------|
| ISO Class 5 | >1 |
| ISO Class 7 | >10 |
| ISO Class 8 or worse | >100 |

If air or surface microbial contamination action levels are reached, taking immediate action will help to quickly eradicate threats and mitigate risks to patient health.



Regulations

Suggested frequency and microbial recovery rates for aseptic processing areas

The USP <1116> - The United States Pharmacopeial Document

Table 2: Suggested Frequency of Sampling for Aseptic Processing Areas

| Clean Room/RABS | |
|---|-----------------------------|
| Sampling Area/Location | Frequency of Sampling |
| Clean Room/RABS | |
| <i>Critical zone (ISO 5 or better)</i> | |
| Active air sampling | Each operational shift |
| Surface monitoring | At the end of the operation |
| <i>Aseptic area adjacent critical zone</i> | |
| All sampling | Each operating shift |
| <i>Other nonadjacent aseptic areas</i> | |
| All sampling | Once per day |
| Isolators | |
| <i>Critical zone (ISO 5 or better)</i> | |
| Active Air Sampling: | Once a day |
| Surface Monitoring: | At the end of the campaign |
| <i>Non-aseptic areas surrounding the isolator</i> | |
| All sampling | Once per month |

The USP <1116> - The United States Pharmacopeial Document

Table 3: Suggested Initial Contamination Recovery Rates in Aseptic Environments

| Room Classification | Active Air Sample (%) | Settle Plate (9cm) 4H Exposure (%) | Contact Plate or Swab (%) | Glove or Garment (%) |
|--|-----------------------|------------------------------------|---------------------------|----------------------|
| Isolator/closed RABS (ISO 5 or better) | <0.1 | <0.1 | <0.1 | <0.1 |
| ISO 5 | <1 | <1 | <1 | <1 |
| ISO 6 | <3 | <3 | <3 | <3 |
| ISO 7 | <5 | <5 | <5 | <5 |
| ISO 8 | <10 | <10 | <10 | <10 |

Ideal for use in less critical environments



Each kit includes:

- Instrument
- Battery charger and cable
- Aspirating head
- Cover head
- LIGHT carrying case
- Calibration certificate



Features:

- Aspirating head with quick bayonet closure system for easy manipulation
- Light weight: Significantly lighter than industry counterparts
- Rechargeable via external port using standard wall charger
- Note: Sampling data must be manually recorded with the MINI

Cat. no.

| | |
|--|---------|
| MINI Kit, 100 liters/min., contact plate | BAS152K |
| MINI Kit, 100 liters/min., Petri plate | BAS153K |
| MINI Kit, 200 liters/min., contact plate | BAS162K |
| MINI Kit, 200 liters/min., Petri plate | BAS163K |

MINI for CompactDry™



TRIO.BAS MINI
designed specifically
for CompactDry™ media.



Each kit includes:

- TRIO.BAS MINI air sampler
- Charging cable
- Stainless steel aspirating head
- Protective cover head
- Carrying case
- CompactDry™ plates sold separately.
See page 35.

Cat. no.

CompactDry Kit, 100 liters/minute

BAS255K

MONO Induction



Ideal for cleanroom use



Each kit includes:

- Instrument
- Induction battery charger and cable
- Aspirating head
- Cover head
- ROBUSTUS carrying case
- Calibration certificate
- IQ, OQ, PQ fillable document

Features:

- No plugs or external connections
- Delayed interval sampling and remote start capabilities
- Data transfer capability via Bluetooth
(Separate purchase required. See page 26)
- **MONO certified to ATEX standards also available**

Cat. no.

| | |
|--|---------|
| MONO Kit, 100 liters/min., contact plate | BAS200K |
| MONO Kit, 100 liters/min., Petri plate | BAS201K |
| MONO Kit, 200 liters/min., contact plate | BAS205K |
| MONO Kit, 200 liters/min., Petri plate | BAS206K |

DUO Induction

Achieve faster sampling times and better statistical results with the DUO!

Ideal for cleanroom use



Each kit includes:

- Instrument
- Induction battery charger and cable
- Aspirating heads
- Cover heads
- ROBUSTUS carrying case
- Calibration certificate
- IQ, OQ, PQ fillable document

Features:

- No plugs or external connections
- Two aspirating heads for sampling on separate plates to determine bacterial and yeast/mold counts
- Programmable for simultaneous or interval aspirating times
- Data transfer capability via Bluetooth. (Separate purchase required. **See page 26**)
- **Ask about DUO certified to ATEX standards**

Cat. no.

DUO Kit, 100 liters/min., contact plate

BAS220K

DUO Kit, 100 liters/min., Petri plate

BAS221K

DUO Kit, 200 liters/min., contact plate

BAS225K

DUO Kit, 200 liters/min., Petri plate

BAS226K

MONO & DUO Cable

WITH SAMPLING DATA TRANSFER CABLE

A solution to capturing data for facilities concerned with data transmissions, or that do not permit wireless transfer of information. Data from the instrument memory can be transferred to a PC via the included transfer cable and dedicated software (see Cat. no. BAS295 or BAS296, required for data download to another device).

TRIO.BAS MONO with Sampling Data Transfer Cable

Each kit includes:

- Instrument
- Stainless steel aspirating head
- Protective cover head
- Charging and Transfer cables
- ROBUSTUS carrying case
- IQ, OQ, PQ fillable document



Cat. no.

MONO Kit, 100 liters/min., contact plate

BAS211K

MONO Kit, 100 liters/min., Petri plate

BAS212K

MONO Kit, 200 liters/min., contact plate

BAS213K

MONO Kit, 200 liters/min., Petri plate

BAS214K

TRIO.BAS DUO with Sampling Data Transfer Cable

Each kit includes:

- Air sampler
- Two stainless steel aspirating heads
- Two protective cover heads
- Charging and Transfer cables
- ROBUSTUS carrying case
- IQ, OQ, PQ fillable document



Cat. no.

DUO Kit, 100 liters/min., contact plate

BAS222K

DUO Kit, 100 liters/min., Petri plate

BAS223K

DUO Kit, 200 liters/min., contact plate

BAS231K

DUO Kit, 200 liters/min., Petri plate

BAS232K

AIRBIO ONE and DUO

The vertical design of AIRBIO instruments is ideal for work surface stability, and monitoring laminar air flow systems.



The AIRBIO DUO holds two culture media plates for sampling with separate media at one time, or to program sampling before (at rest), during (in operation), and at the end of each processing cycle.

Each kit includes:

- Air sampler
- Bluetooth capability
- Stainless steel aspirating head(s)
- Stainless steel protective aspirating head cover(s)
- Charging cable
- Data transfer cable
- Robustus carrying case
- Factory calibration certificate
- IQ OQ PQ fillable documents

Features:

- Compliant according to EN/ISO 14698-1, GMP, GLP, USP, and data integrity CFR21 and GAMP-5
- Volume of aspirated air is either 100 liters per minute or 200 liters per minute
- Sample air intake volume programs from 30-2000 liters with 17 preset programs
- Use either standard culture media plates: Petri plates or contact plates
- IP65 rated protection from dust and water



Cat. no.

| | |
|--|---------|
| AIRBIO MONO 100 liters/min contact plate | BAS445K |
| AIRBIO MONO 100 liters/min Petri plate | BAS446K |
| AIRBIO MONO 200 liters/min contact plate | BAS447K |
| AIRBIO MONO 200 liters/min Petri plate | BAS448K |
| AIRBIO DUO 100 liters/min contact plate | BAS479K |
| AIRBIO DUO 100 liters/min Petri plate | BAS480K |
| AIRBIO DUO 200 liters/min contact plate | BAS481K |
| AIRBIO DUO 200 liters/min Petri plate | BAS482K |

AIRBIO ONE RAPID-VIRUS

*For airborne viable particle sampling.
It is specifically designed for total
pathogen surveillance of bacteria,
fungi, yeast and viruses.*



Options:

- Use to collect samples in liquid for subsequent rapid analytical identification steps by PCR, etc.
- Use with the traditional impact on agar culture media plate method to count the colonies (CFU).

Features:

- Collection unit is completely sterilizable
- The tri-clamp stainless steel tube system facilitates cleaning and sterilization
- Liquid sampling permits rapid methods (e.g.:PCR), while traditional agar sampling requires
- incubation of media plates to obtain sample results.

Cat. no.

AIRBIO ONE RAPID-VIRUS instrument

BAS2448K

Workflow:



1. Remove the protective cover head from AIRBIO



2. Connect the sterile virus system to the aspirating chamber of AIRBIO



3. Disconnect the protective bottle



4. Connect the bottle with the collection liquid



5. The air sampler is ready



6. Start sampling



7. The sample is transferred to the laboratory



8. The sample is analysed via PCR System

Is a HEPA filter-equipped viable air sampler necessary for monitoring cleanrooms?

Active microbial air samplers for viable contamination monitoring in cleanroom isolators and RABS (Restricted Access Barrier Systems) should be designed to avoid contamination risk from the instrument.

The electrical motor of TRIO.BAS™ instruments does not produce particles if the manufacturer recommended cleaning procedures are followed during sampling preparation. Furthermore, expelled air flows out from the aspirating chamber without any contact with the internal parts of the instrument.

When permanently located in cleanroom isolators and RABS, TRIO.BAS™ instruments without HEPA filters are used without concern of contamination from particles.



Cat. no.

TRIO.BAS MONO
Filter Kit, 100 liters/min.,
contact plate

BAS171K

TRIO.BAS MONO
Filter Kit, 100 liters/min.,
petri plate

BAS170K



Cat. no.

AIRBIO ONE Filter Kit,
100 liters/min.,
contact plate

BAS450K

AIRBIO ONE Filter Kit,
100 liters/min., petri plate

BAS449K

If the viable air monitoring plan involves transferring the air sampling instrument from different ISO classification environments with varied levels of microbial contamination (ISO Class 7 to ISO Class 5), a risk of cross-contamination is possible.

Under such conditions, a HEPA filter equipped TRIO.BAS™ instrument, or a HEPA filter satellite option for use with RABS/ISOLATOR instrument formats, is suggested to avoid the transfer of particles from the inside of the aspirating chamber.

TRIO.BAS™ HEPA filter formats include a replaceable HEPA filter fixed to the air outlet of the instrument or satellite to stop the potential transfer of particles.

Satellite for use with BAS268K, BAS269K, BAS270K, BAS271K



Cat. no.

Stainless steel satellite
with HEPA FILTER,
contact plate

BAS262K

Stainless steel satellite
with HEPA FILTER,
Petri plate

BAS263K



A smoke/exhaust study confirms
the expelled air is not re-aspirated
by the same air sampler.

Cat. no.

Replacement HEPA filters,
10pk

BAS531

RABS ISOLATOR and MULTIFLEX Systems

Specifically designed to meet Pharmaceutical and Biopharmaceutical Cleanroom Regulatory Requirements and Standards



**Data transfer requires additional purchase of BAS296 software*

RABS ISOLATOR

External Command Unit and One Satellite

with the ability to add up to two additional satellite kits

Each kit includes:

- Instrument
- Satellite aspirating head
- Stainless steel head cover
- Five meter long satellite attachment cable
- Charging and Data transfer* cable
- Robust carrying case
- Calibration certificate
- IQ OQ PQ fillable guideline documents

Cat. no.

| | |
|--|---------|
| RABS ISOLATOR 1 Satellite Pack, 100 liters/min., contact plate | BAS268K |
| RABS ISOLATOR 1 Satellite Pack, 100 liters/min., Petri plate | BAS269K |
| RABS ISOLATOR 1 Satellite Pack, 200 liters/min., contact plate | BAS270K |
| RABS ISOLATOR 1 Satellite Pack, 200 liters/min., Petri plate | BAS271K |

Standard Satellite Kits with Cables

Each kit includes:

- Satellite air sampling chamber
- Stainless steel aspirating head
- Stainless steel head cover
- Five meter long satellite attachment cable
- LIGHT carrying case

| | |
|---------|---------|
| Contact | BAS260K |
| Petri | BAS261K |



HEPA-Equipped Satellites

- Add one to three satellites
- Stainless steel HEPA satellite
- Stainless steel aspirating head
- Stainless steel head cover
- Five meter long satellite attachment cable
- One box of HEPA filters

| | |
|---------|---------|
| Contact | BAS262K |
| Petri | BAS263K |




Cannot be used concurrently with standard, or other, satellite configurations. Not for use with MULTIFLEX.

Satellites with Tri-Clamps

- Satellite and Tri-Clamp
- Aspirating heads and covers sold separately (BAS330/331 and BAS465)

| | |
|---------|---------|
| Contact | BAS320K |
| Petri | BAS321K |

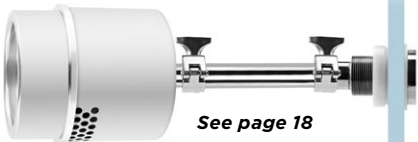

See page 18

- Compliant according to EN/ISO 14698-1, GMP and GLP
- Built in ISO 9001 Premises
- Data Integrity per CFR 21 and GAMP5
- IQ OQ PQ Guidance Documents Included
- Fabricated in AISI 316 Stainless Steel
- Multiple satellites allow use of separate media at the same time to determine bacterial and fungal counts
- Bluetooth data transfer capability
- Monitor separate cleanrooms with a single external command unit and satellites attached to 5 meter long cables
- Satellite designs allow sampling inside RABS or Isolators
- Optional 20 meter long satellite cables available (BAS272)
- Sample up to 70,000 liters on a single battery charge



**Data transfer requires additional purchase of BAS296 software*

MULTIFLEX 1

Command Unit + Built-in Aspirating Head
with the ability to add up to two additional standard satellite kits

Each kit includes:

- Instrument
- Stainless steel head cover
- Charging and Data transfer* cable
- Robustus carrying case
- Calibration certificate
- IQ OQ PQ fillable guideline documents

Cat. no.

| | |
|---|---------|
| MULTIFLEX 1, 100 liters/min., contact plate | BAS483K |
| MULTIFLEX 1, 100 liters/min., Petri plate | BAS484K |
| MULTIFLEX 1, 200 liters/min., contact plate | BAS485K |
| MULTIFLEX 1, 200 liters/min., Petri plate | BAS486K |



**Data transfer requires additional purchase of BAS296 software*

MULTIFLEX 1+2

Command Unit + Built-in Aspirating Head and Two Independent Satellites


Each kit includes:

- Instrument
- Three stainless steel aspirating heads with stainless steel head covers
- Two 5 meter long satellite attachment cables
- Charging and Data transfer* cable
- Calibration certificate
- IQ OQ PQ fillable guideline documents

Cat. no.

| | |
|---|---------|
| MULTIFLEX 1+2, 100 liters/min., contact plate | BAS474K |
| MULTIFLEX 1+2, 100 liters/min., Petri plate | BAS475K |
| MULTIFLEX 1+2, 200 liters/min., contact plate | BAS476K |
| MULTIFLEX 1+2, 200 liters/min., Petri plate | BAS477K |

Remote Aspirating Tube System

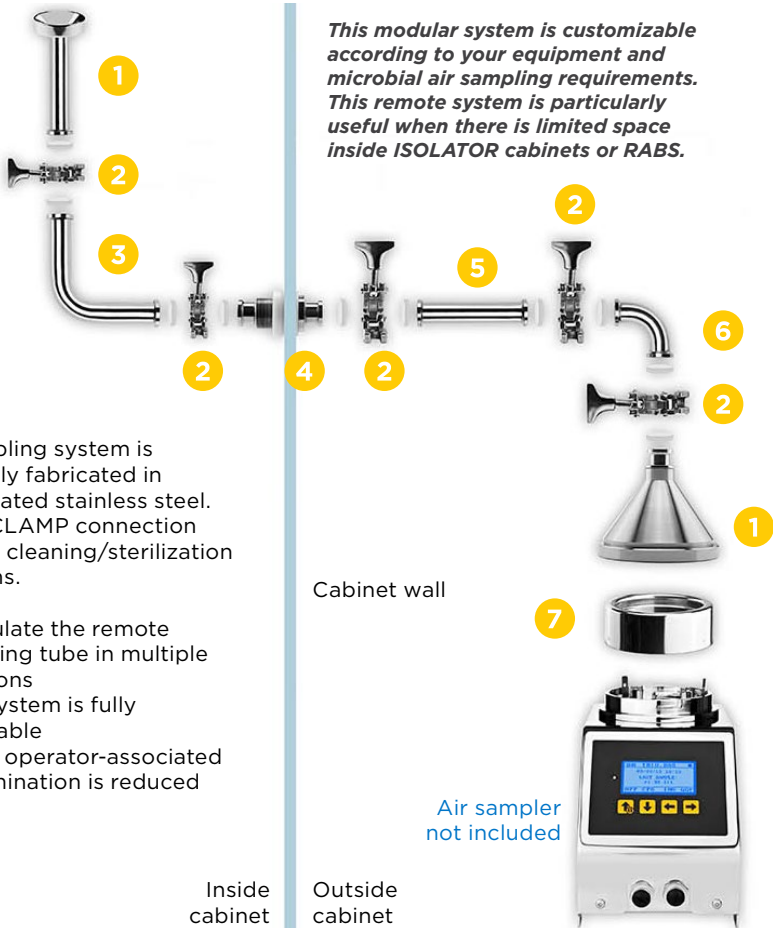


This modular system is customizable according to your equipment and microbial air sampling requirements. This remote system is particularly useful when there is limited space inside ISOLATOR cabinets or RABS.

This sampling system is completely fabricated in AISI 316 rated stainless steel. The TRI CLAMP connection simplifies cleaning/sterilization operations.

- Manipulate the remote aspirating tube in multiple directions
- Tube system is fully sterilizable
- Risk of operator-associated contamination is reduced

| | | Cat. no. |
|---|---|----------|
| 1 | REMOTE TUBE SYSTEM, stainless steel aspirating tube and end clamp with stainless steel aspirating bell and end clamp to connect the system to Air Sampler | BAS181 |
| 2 | Stainless steel tri-clamp with silicone gasket | BAS193 |
| 3 | Stainless steel 90° tube with 2 clamp ends, size wheelbase 88, 90mm, radius 25, 40mm | BAS195 |
| 4 | Stainless steel wall connection with 2 clamp ends | BAS198 |
| 5 | Stainless steel tube with 2 clamp ends, 10cm | BAS194 |
| 6 | Stainless steel 90° tube with 2 clamp ends, size of wheelbase 40mm | BAS188 |
| 7 | Stainless Steel Aspirating Head, Petri plate, shown as part of the MULTIFLEX 1 microbial air sampler | BAS331 |

Remote Aspirating Head System

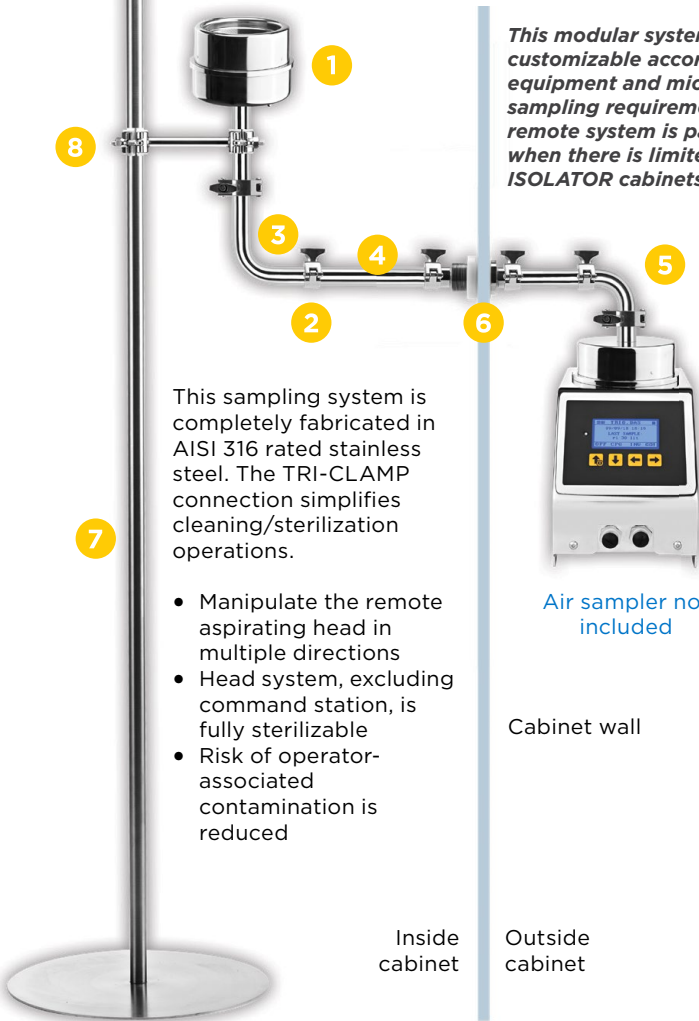

This modular system is customizable according to your equipment and microbial air sampling requirements. This remote system is particularly useful when there is limited space inside ISOLATOR cabinets or RABS.

This sampling system is completely fabricated in AISI 316 rated stainless steel. The TRI-CLAMP connection simplifies cleaning/sterilization operations.

- Manipulate the remote aspirating head in multiple directions
- Head system, excluding command station, is fully sterilizable
- Risk of operator-associated contamination is reduced

Air sampler not included

Cabinet wall

Inside cabinet

Outside cabinet


Cat. no.

| | | |
|---|--|------------|
| 1 | Stainless Steel Contact or Petri format Aspirating Head | BAS191/192 |
| 2 | Stainless Steel tri-clamp with silicone gasket | BAS193 |
| 3 | Stainless Steel elbow 90° tube with 2 clamp ends - size wheelbase 88,90 mm, radius 25,40mm | BAS195 |
| 4 | Stainless Steel tube with 2 clamp ends - 10cm length | BAS194 |
| 5 | Stainless Steel elbow short tube 90° with 2 clamp ends - size wheelbase 40mm | BAS188 |
| 6 | Stainless Steel wall connection with 2 clamp ends | BAS198 |
| 7 | Stainless Steel floor pole - diameter 25 cm, 1 mt height | BAS199 |
| 8 | Stainless Steel hexagonal pipe hook with gaskets and adapter | BAS197 |

TRI CLAMP Satellite

A customizable air sampling satellite option for use with RABS ISOLATOR or MULTIFLEX instruments. Ideal for use in cleanroom barrier systems where space is limited.

- For integration with the RABS Isolator or MULTIFLEX instrument
- Compact design for use in RABS or isolators
- Simple electric cable
- Fabricated in AISI 316 rated stainless steel
- Easily mounted horizontally or vertically
- Tri-clamp coupling system for easy attachment and disassembly during maintenance and calibration
- Electric connection cable to the control unit remains protected inside the stainless steel tube
- Use either Petri or contact plates. Adapts to ready-to-use technopolymer Daily Shift Heads

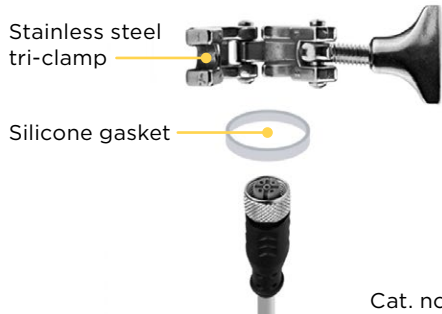

Cat. no.

| | |
|--|---------|
| Stainless steel tri-clamp aspirating head satellite for Contact plate, with stainless steel head cover and cable connector | BAS320K |
| Stainless steel tri-clamp aspirating head satellite for Petri plate, with stainless steel head cover and cable connector | BAS321K |
| Stainless steel tri-clamp with silicon gasket | BAS193 |
| Stainless steel tube with 2 clamp ends, 10cm length | BAS194 |
| Stainless steel wall connection with retaining nut, for walls max of 5mm thick | BAS198 |
| Stainless steel wall connection with retaining nut, for walls over 5mm thick | BAS189 |



Designed to test the presence of microorganisms in compressed air or gas lines.

- Calibrated regulator guarantees 100 liters per minute flow rate
- Autoclavable with no disassembly required
- Automated to end sample after programmed volume is reached*
- Signal lets operator know sample is complete*
- Sampling data is retrievable for download via proprietary printer or dedicated software (MONO) or to record manually (MINI).
- No glass valves or meters to crack or break

*When used with MONO or MINI instrument



Each kit includes:

- Stainless steel electro-valve
- Gas connection
- Stainless steel fixing system for air sampler
- MINI or MONO unit
- Induction battery charger and cable
- Aspirating head
- Head cover
- Two IQ, OQ, PQ fillable documents
- LIGHT carrying case

Cat. no.

| | |
|--|---------|
| TRIO.GAS System + MONO air sampler, 100 liters/min., contact plate | BAS650K |
| TRIO.GAS System + MONO air sampler, 100 liters/min., Petri plate | BAS651K |
| TRIO.GAS System + MINI air sampler, 100 liters/min., contact plate | BAS654K |
| TRIO.GAS System + MINI air sampler, 100 liters/min., Petri plate | BAS655K |


TRIO.GAS System Aspirating Chamber Kit

Each kit includes:

- Stainless steel electro-valve
- Gas connection
- Stainless steel fixing system for air sampler
- Aspirating head
- Cover head
- Carrying case
- IQ, OQ, PQ fillable documents
- MONO or MINI NOT included

Important Note:

If manually timed samples and hand recorded sampling data is acceptable, the MINI or MONO instrument is not necessary. Purchase the GAS system BAS652K or BAS653K, which includes aspirating gas chambers. Choose kit based on preferred plate size.

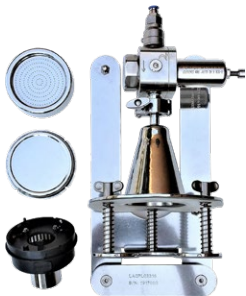


Cat. no.

| | |
|---|---------|
| TRIO.GAS System Aspirating chamber kit, contact plate | BAS652K |
| TRIO.GAS System Aspirating chamber kit, Petri plate | BAS653K |

FLO.GAS

Test for the presence of microorganisms in compressed air and various gases supplied from tanks and pipes under pressure. The flow rate regulator valve must be adjusted prior to sampling.

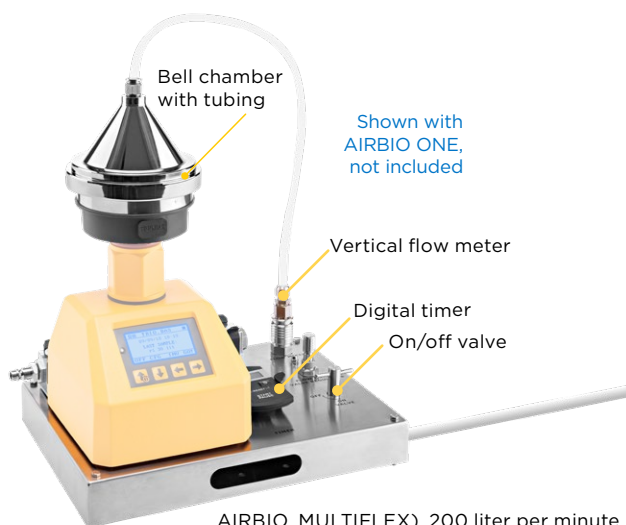

- Suitable for testing compressed air, nitrogen, CO₂ and Argon (Use a separate flow meter, BAS596, for CO₂ and/or Argon testing)
- Autoclavable
- Compatible with ASPI GAS CHAMBER (included), or 100 liter per minute air intake rate air samplers (TRIO.BAS,

AIRBIO, MULTIFLEX). 200 liter per minute samples require a separate air flow meter, Cat. no. BAS595

- Annual calibration recommended

Cat. no.

FLO.GAS System, with flow meter, ASPIGAS chamber, Petri plate aspirating head and digital timer


BAS597

FLO.GAS System, with flow meter, ASPIGAS chamber, contact plate aspirating head and digital timer

BAS598

VERI.GAS

Two performances in one instrument: Test for the presence of microorganisms in compressed air/gas and verify the calibrated air flow rate is within specification. The flow rate regulator valve must be adjusted prior to sampling.

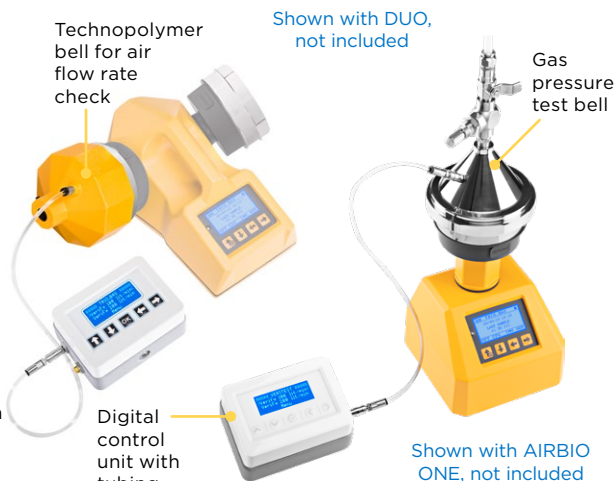

- Suitable for testing compressed air and nitrogen
- Autoclavable
- Compatible with 100 liter per minute air sampler instruments (TRIO.BAS, AIRBIO, MULTIFLEX)
- For gas testing, use the AISI rated stainless steel bell chamber, valve, and regulator

- For calibration verification, use the technopolymer bell chamber
- Annual calibration recommended

Cat. no.

VERI.GAS, for compressed gas test and air flow rate check - with digital control unit, s/s bell chamber for gas testing, technopolymer bell chamber for calibration specification check, connection tube and robust carrying case

BAS599

Calibration Control Equipment

Proper calibration of air sampling devices is a requirement per cGMP and the USP. Hardy Diagnostics' ISO 17025 certified lab is the factory trained and sole North American calibration service center for TRIO.BAS air sampling instruments.

For detailed information, please see the **TRIO.BAS Instrument Services** tab located under **Product Support** on our website:
HardyDiagnostics.com/air-sampler-calibration/

SELF TEST

A manual verification instrument to check the precision of air flow rate.

For control purposes only; not to take the place of annual, official calibrations. The base station induction battery charger with user **SelfTest** is equipped with a system that allows, regardless of auto-calibration already present in the instrument, the ability to check the state of precision of the air flow.


Only compatible with induction battery charging system instrument formats.

Calibration System includes:

- Calibration bell
- Connector cable and charger with outlet plug

For 100 liter/min.
flow instrument BAS351

For 200 liter/min.
flow instrument BAS352

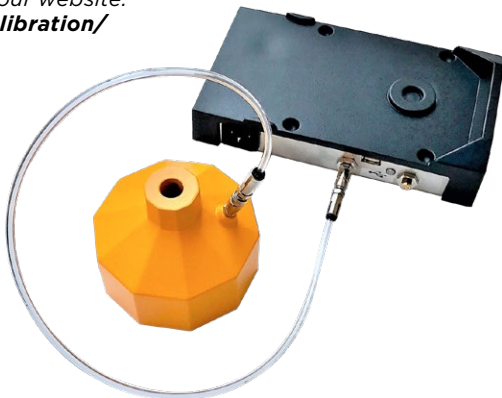

VERITEST

A manual verification instrument to check the precision of air flow rate.

For control purposes only; not to take the place of annual, official calibrations. Suitable for all samplers and satellites.

Each BAS353



DAILY SHIFT HEADS

Individually packaged, sterile by irradiation, and includes a certificate of sterility assurance.

For use on the MINI, MONO, DUO, TRIO and MULTISTATION units.

Technopolymer aspirating heads are ready for immediate cleanroom use.



Benefits of Daily Shift Heads:

- A daily certificate of sterilization assurance is requested by regulatory authorities as part of the usual auditing process. Use of DAILY SHIFT aspirating heads eliminates the cost of sterilization after each sample cycle and the consequent task of preparing the certificate of sterility
- DAILY SHIFT heads offer convenience and time savings during periods of unusually heavy workloads
- DAILY SHIFT heads offer peace of mind in case of autoclave breakage or unavailability
- DAILY SHIFT heads can be used on the same day, on the same shift, inside the same cleanroom
- Semi-transparent to confirm agar plate has been properly inserted
- Six year shelf life from date of sterilization

Cat. no.

DAILY SHIFT HEAD, contact plate, 27/bx

BAS340

DAILY SHIFT HEAD, Petri plate, 27/bx

BAS341

Data Transfer Options

For TRIO.BAS instruments with Bluetooth or cable Data Transfer Capability, choose the preferred data download method. A separate purchase is required.

Tablet Device*

(for data transfer)

Portable tablet for remote control of instrument(s), data download storage and transfer capabilities.

**Requires software BAS296*

Each

BAS301



Bluetooth Printer

Small footprint. Pair with instrument for immediate printout of sampling data, or choose to store data and print out sampling history at a later date.

Bluetooth printer, 11 x 8 x 5

BAS520

Paper rolls (57mm, 10/bx)

BAS421



Bluetooth Key*

For transferring sampling data to a PC that does not include built-in Bluetooth feature.

**Requires software BAS296*

Each

BAS420



Data Transfer Software*

Biological Air Sampler Software (BAS)**

*One software license per PC is required. Multiple instruments can be used with a single software license.

Designed to facilitate paperless record keeping of all aspects of environmental monitoring, BAS software is the ideal TRIO.BAS instrument accompaniment for cleanroom environments.

This innovative software tool facilitates the collection and analysis of air sampling data, as well as tracability of the complete environmental monitoring sampling process. Include pictures of the CFU count from sampled media plates when the CFU camera instrument, BAS337 (pg 29), is incorporated.

- Make the transition from paper to paperless records and enhance standardization with BAS Software
- Simplifies the entry of user profiles and sampling schedules

Each BAS296

***For Bluetooth capable Windows systems*



Data Management Air Sampler Software (AS)**

AS Software is a simple data transfer solution, suggested when the main objective is to transfer data to a PC for record keeping or analysis. AS software does not meet Data Integrity or ALCOA regulatory standards.



***For Bluetooth capable Windows systems*


- Transfer data from Bluetooth or data transfer cable-equipped instrument, or from the proprietary TRIO.BAS tablet, to a PC via Bluetooth

Each BAS295

ALCOA Compliant Records of Passive, Active, Compressed Air, & Surface Samples (P.A.C.A.S.) are possible with BAS296 Software

BAS Software for managing active and passive air, compressed gas and surface sampling cycles

Final report including plate pictures (air, compressed gas and surface) before and after colony count

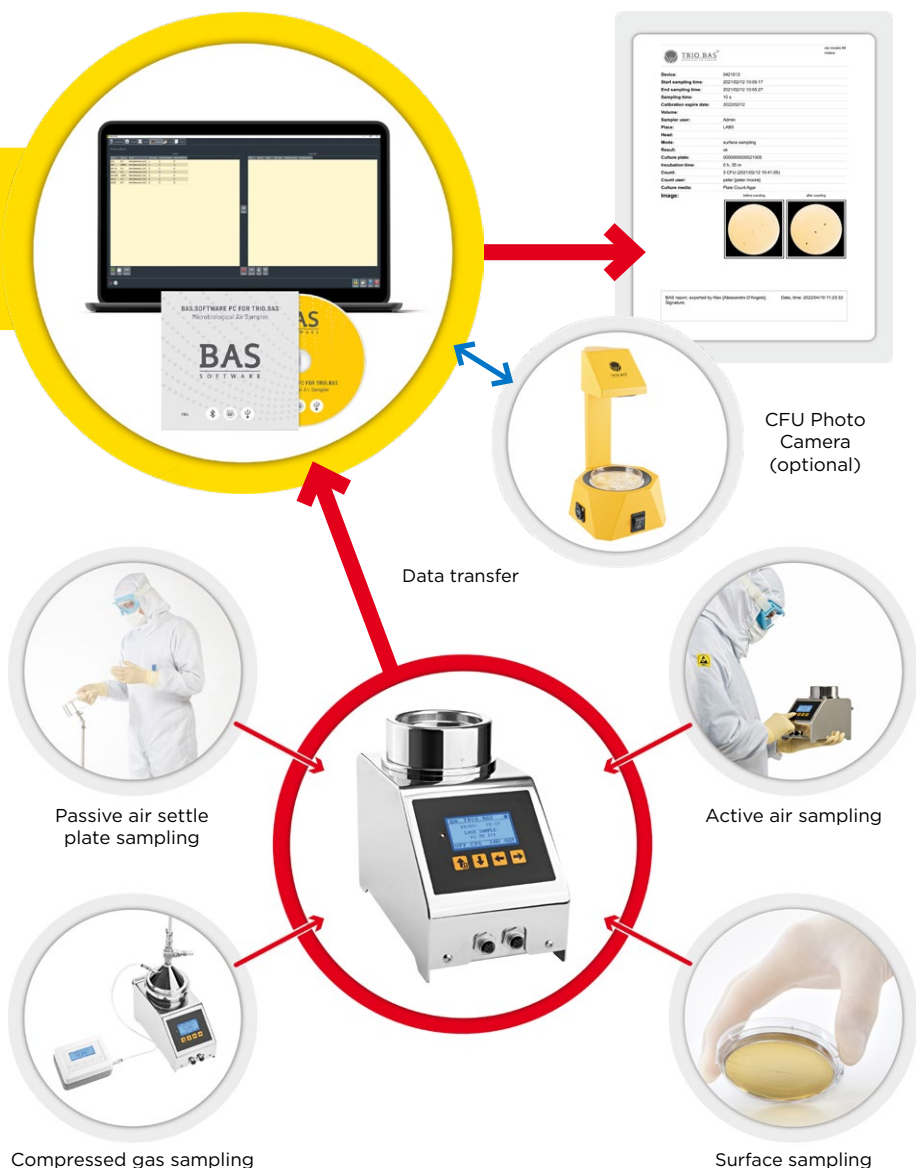

Comparing AS SOFTWARE and BAS SOFTWARE



| COMPARISONS | AS SOFTWARE BAS295 | BAS SOFTWARE BAS296 |
|--|------------------------------------|--|
| Password protected login | X | ✓ |
| Connectivity of TRIO.BAS to PC via Bluetooth/cable | ✓ | ✓ |
| Remote control of TRIO.BAS instrument | ✓ | ✓ |
| Protected password report exporting | Free export file .pdf .csv .asd | Password export file .pdf .csv .bas .xml |
| Configuration of users and places (sample sites) via software | X | ✓ |
| Track and report on complete sampling process (P.A.C.A.S.) | X | ✓ |
| Include photos of sampled agar plates in reports when utilizing the CFU photo camera | X | ✓ |
| Compliant to FDA 21 CFR part II data integrity and ALCOA regulatory standards | X | ✓ |
| Audit trail | X | ✓ |

CFU Photo Camera

Obtain culture media plate photos immediately before and after the manual CFU count with the CFU Photo Camera. The CFU Photo Camera integrates with BAS296 software to capture unalterable records, per 21 CFR part 11, of air sampling media plate count results, with the ability to export data to a pdf file or printed paper copy.



- Light System: LED technology
 - Camera resolution: full HD 1080p
 - Construction: technopolymer antibacterial treatment
 - Easy to connect to the PC
 - Works with BAS296 SOFTWARE only
- Each BAS337



Before counting



After counting

Accessories

Technopolymer Standup Holder

for MINI

Each BAS370

for MONO, DUO
and TRIO

Each BAS376

for AIRBIO DUO

Each BAS377



Support stand for MULTIFLEX

A stainless steel support stand for stability and safety when positioning MULTIFLEX instruments to best capture laminar flow.

Each BAS379



Stainless Steel Vertical Hook and Knob

For table, cart, or tripod. Adaptable to MINI, MONO, DUO, and TRIO.

12x15x25 cm.

Each BAS521

Stainless Steel Wall/ Table Holder

For wall mount or table.

Each BAS530



Stainless Steel AISI 316 Tripod

150cm to 200cm. Ideal for cleanroom use. Device not included.

Each BAS523



Maxi Floor Tripod

The sturdy MAXI Tripod goes to great lengths to facilitate sampling at heights up to 366cm (12 ft).

Each BAS387

Center Ball Head

An optional head for MAXI Tripod the to rotate the air sampler in different directions. Not recommended when extending the MAXI Tripod to its maximum height.

Each BAS388



Soft Carrying Case

Soft carrying case with tote handles for transporting the TRIO.BAS™ MINI or MONO.

Each BAS390



Multi Holder Carts


Mini Multi Holder Cart

Stainless steel AISI 316 cart on wheels equipped with three shelves. Ideal inside cleanroom.

- Dimensions: 25x35x70cm


Each

BAS371

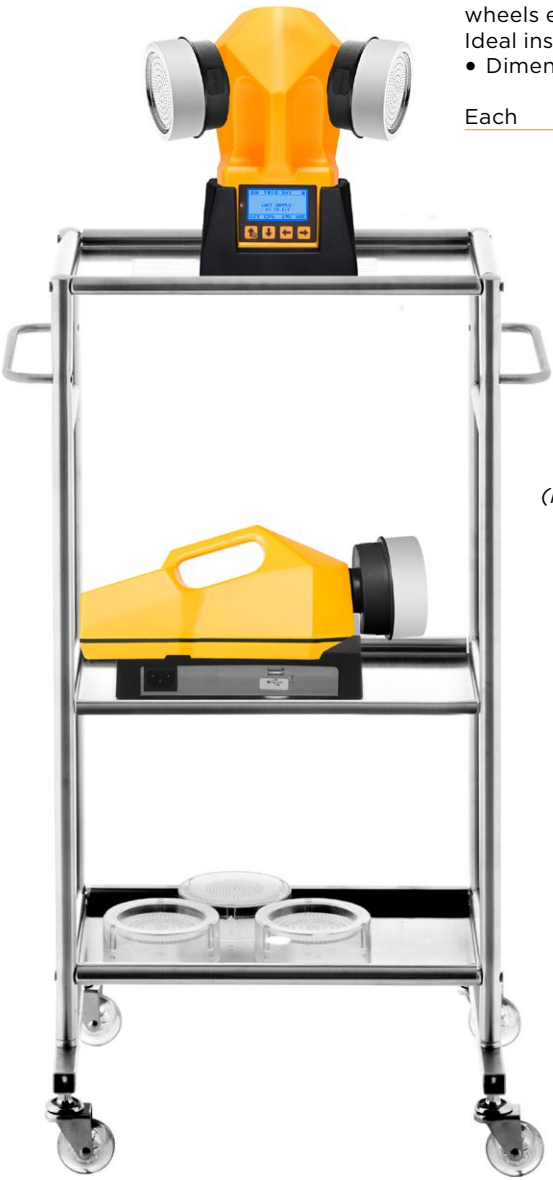

(Pictured with TRIO, MONO and stand, not included)

Multi Holder Carts

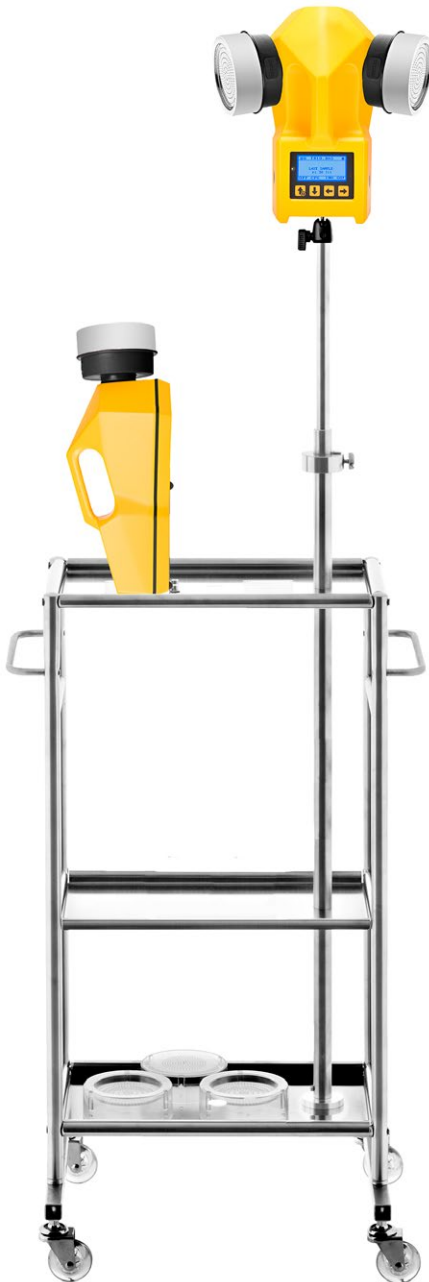

Maxi Multi Holder Cart

Stainless steel AISI 316 cart on wheels with fixed vertical hook holder and adjustable vertical extension post. Ideal inside cleanroom.

- Dimensions: 25x35x70cm
- Extension post adjustable height: 100cm to 210cm

Each


BAS372

(Pictured with DUO and MONO, not included)

TRIO.SETTLE

Ideal for settling
plate exposure under
unidirectional air flow

Standardize your passive air monitoring protocols with TRIO. SETTLE. Fully autoclavable AISI 316 rated stainless steel support for Settling Plates.

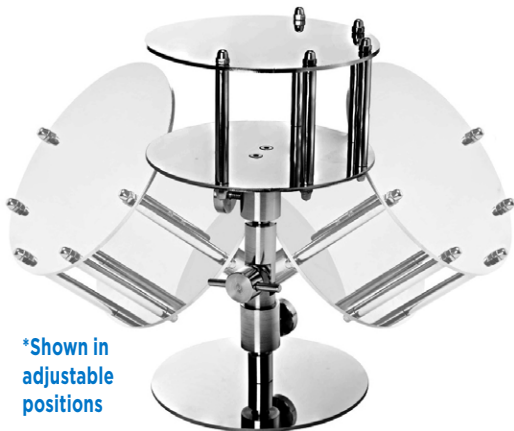

***Shown in adjustable positions**

The TRIO.SETTLE ideally standardizes the position of the culture plate reducing the risk of contamination during sampling.

Position culture media (settling plates) for best articulation to provide a quantitative analysis of airborne microorganisms deposited over a set period of exposure.

TRIO.SETTLE features:

- AISI316 stainless steel fabrication
- Upper surface disc to accommodate an open culture plate during sampling
- Lower surface disc to accommodate the lid of the culture plate
- Adjustable to obtain different inclinations to the agar surface related to the direction of unidirectional airflow, avoiding laminar flow turbulence.
- Floor base model includes floor base and a 1 meter tall stainless steel tripod

Cat. no.

TRIO.SETTLE table plate stand

BAS367

TRIO.SETTLE floor plate stand

BAS368



TRIO.BAS IN-REST Plate Holder

To facilitate operations and reduce contamination risks during sampling activities. Useful as a temporary support for culture media plates and/or lids, as well as TRIO.BAS instrument sampling heads.

Each

BAS180

Cleanroom Bag

Sterile bags ideal for safely containing and transporting Petri dishes, Contact plates, swabs, or other objects, outside of the cleanroom.

5/pk BAS381BX

750 Bags/Case BAS381



Prepared Culture Media and CompactDry™



Irradiated, Triple-Bagged, Mono Plates

For air sampling and settling plates

Tryptic Soy Agar (TSA) with Lecithin and Tween® 80, SterEM™, USP

For use as a general growth medium for establishing microbiological trends, alerts and action levels in biologically controlled environments.

15x100mm plate, 26ml fill, 10/pk

W520

SabDex (Sabouraud Dextrose) Agar, SterEM™, USP

For the cultivation of fungi. Irradiated, 15x100mm plate, triple bagged, 10/pk

W565

Irradiated, Triple-Bagged, Contact Plates

For air and surface sampling

TSA with Lecithin and Tween® 80, USP

For the cultivation and enumeration of microorganisms with Lok-Tight™ lid. Irradiated, triple bagged, 10/pk

P520

SabDex Agar with Lecithin & Tween® 80, USP

For the cultivation of fungi. With Lok-Tight™ lid. Irradiated, triple bagged, 10/pk

P595



CompactDry™ Media

For TRIO.BAS CompactDry Instrument Only.

See page 6

Total Plate Count

CompactDry™ TC, Total Plate Count, 240/box.

54081

Yeast & Mold

CompactDry™ YM, Yeast and Mold, 240/box

54083

Yeast & Mold Rapid

CompactDry™ YMR, Yeast Mold Rapid, 240/box

54084





TRIO.BASTM
BIOLOGICAL AIR SAMPLER



 **HARDY**
DIAGNOSTICS
A Culture of Service™



TRIO.BASTM
BIOLOGICAL AIR SAMPLER



Headquarters

1430 West McCoy Lane
Santa Maria, CA 93455
800.266.2222
Sales@HardyDiagnostics.com
HardyDiagnostics.com

Distribution Centers

Santa Maria, California
Olympia, Washington
Salt Lake City, Utah
Phoenix, Arizona
Dallas, Texas
Des Moines, Iowa
Springboro, Ohio
Lake City, Florida
Albany, New York
Raleigh, North Carolina

Copyright © 2025 Hardy Diagnostics



Hardy Diagnostics
donates 1% of net
profits to charity.



FM 572526

Hardy Diagnostics has a Quality Management System that is certified to ISO 13485 and is a FDA licensed medical device manufacturer.



071825ss