




RAPID AGGLUTINATION

TEST KITS

RAPID AGGLUTINATION TEST KITS

Contents

- | | | | |
|---|-------------------------|---|--------------------------------|
| 3 | ASO (Anti-Streptolysin) | 4 | Rheumatoid Arthritis |
| 3 | <i>Cryptococcus</i> | 5 | Rubella |
| 3 | <i>E. coli</i> O157 | 5 | <i>Staphylococcus</i> and MRSA |
| 4 | Mononucleosis | 6 | Strep Grouping |

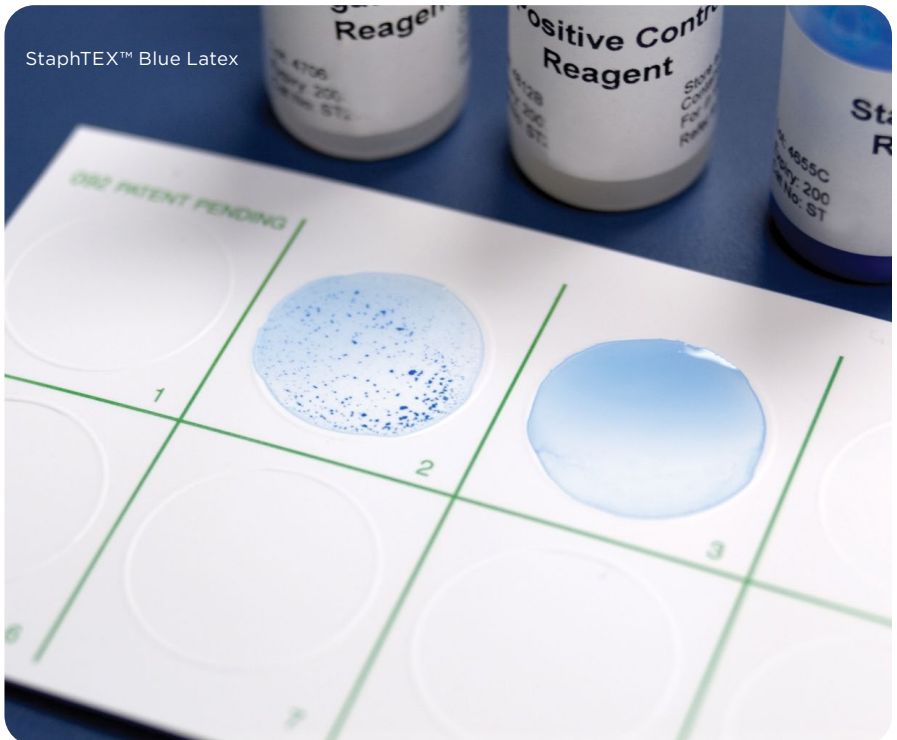


Agglutination Testing

Agglutination testing is a technique used in microbiology and immunology to detect and identify certain antigens or antibodies by observing the clumping (agglutination) of particles. This method is commonly used in the identification of microorganisms, especially bacteria, and in the detection of antibodies in blood serum.

The principle behind agglutination testing is based on the specific interactions between antigens and antibodies. Antigens are molecules that can stimulate an immune response, and antibodies are proteins produced by the immune system in response to the presence of antigens.

In microbiology, agglutination tests are often used to identify and classify bacteria based on the specific antigens they possess.

Hardy Diagnostics is proud to showcase a variety of rapid agglutination test solutions to aid our laboratory partners in diagnosing infectious diseases, and identifying bacterial strains. Our testing solutions have long shelf lives and most provide answers to you in just a few minutes or less.

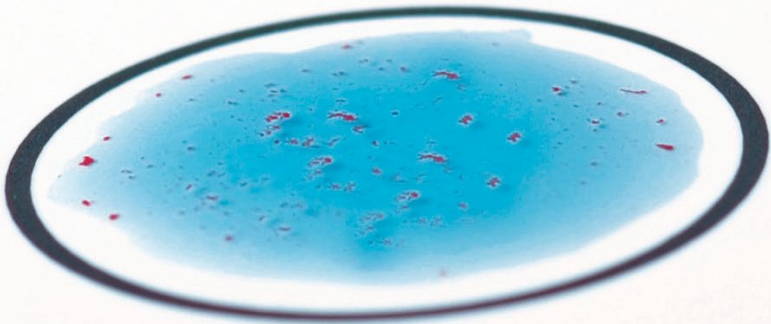

ASO (Antistreptolysin-O)

ASO Slide Test

A latex slide agglutination assay for the qualitative and semiquantitative detection of antistreptolysin-O (ASO) antibodies in human serum. No initial dilution of patient sample is required for this test. Contains reactive and nonreactive controls.

- Sensitivity: 200 IU/ml
- Specificity: >95%
- Results in two minutes

50 tests	300050
100 tests	300100

ASO Control Kit

If needed, additional controls can be purchased seperately.

2/pk	305001
------	--------

CRP Slide Test

A latex slide agglutination assay the qualitative and semiquantitative detection of C-reactive protein (CRP) in human serum. No initial dilution of patient samples is required for this test. Contains reactive and nonreactive controls.

- Sensitivity: 6mg/L
- Specificity: >95%
- Results in two minutes

50 tests	500050
----------	--------

Cryptococcus

Cryptococcus Antigen Test Kit

A qualitative or semi-quantitative, 5-minute latex agglutination test for detection of polysaccharide antigens associated with *Cryptococcus neoformans* infection. This test can be performed from serum or CSF specimens. Contains both positive and negative controls.

- Results in five minutes

50 tests	R30851501
----------	-----------

E. coli O157

E. coliPRO™ O157 Latex

A latex agglutination test for the presumptive detection of *E. coli* serogroup O157 antigen from colonies isolated from culture media. Proprietary antibodies in this test minimize cross-reactions with *Escherichia hermannii* and other sorbitol-negative *Escherichia* spp. Test contains both positive and negative controls.

- Results in two minutes

50 tests	PL070HD
----------	---------

E. coliPRO™
O157 Latex



Mononucleosis

Color Slide II Mononucleosis Test


A color-enhanced slide agglutination test for detection of heterophile antibodies associated with infectious mononucleosis in serum, plasma, or finger-tip (whole) blood. Test is both qualitative and semiquantitative. Contains both positive and negative controls.

- Results in one minute

50 tests

2468944

Color Slide II
Mononucleosis Test

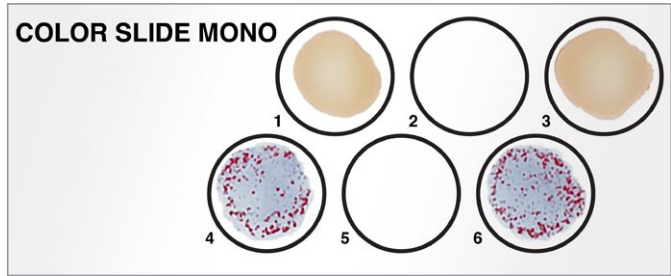

Color Mono II Test

An agglutination test for the qualitative and semi-quantitative detection of heterophile antibodies, in human serum or plasma, associated with infectious mononucleosis. No initial dilution of patient samples is required for this test; reactive and nonreactive controls included.

- 18-24 month refrigerated shelf life from date of manufacture
- Sensitivity: 1:28 Davidsohn titer
- Specificity: >95%
- Results in two minutes

50 tests

450050

Rheumatoid Arthritis

RF Direct Slide Test

The ASI RF Direct Slide Test is a slide agglutination assay for the qualitative and semiquantitative determination of rheumatoid factor (RF) in human serum. No initial dilution of patient samples is required for this test; reactive and nonreactive controls included.

- Results in two minutes

100 tests

700100DC



Color Mono II Test

Rubella

Rubella Latex

A rapid latex agglutination test for the qualitative and semi-quantitative determination of rubella virus antibodies in serum. This test aids in the diagnosis of recent or active rubella infection and in the determination of immune status. Results from this test correlate with results from hemagglutination inhibition assays. Reactive and nonreactive controls are included.

- Results in eight minutes
- Sensitivity: 1-2 IU/ml (undiluted); 10-20 IU/ml (diluted)

100 tests 600100

Staphylococcus and MRSA

StaphTEX™ Blue Latex

A rapid latex agglutination slide-card test for the identification of *Staphylococcus aureus* from culture media. Tests for both coagulase and protein A found in *S. aureus*. Coagulase is detected in both free and bound form (clumping factor). Contains positive and negative controls.

- Results in less than 30 seconds
- Blue latex on white card for easy read-out

50 tests ST50

200 tests ST200

1,000 tests ST1000

MASTALEX™ MRSA

A rapid slide latex test for the detection of penicillin-binding protein 2' (PBP2' also called PBP2a) in isolates of *Staphylococcus*. This test assists in identifying of methicillin-resistant *Staphylococcus aureus*. Contains negative controls.

- Results after three minutes
- May be used for screening or diagnosis

48 tests RST501



StaphTEX™
Blue Latex

Strep Grouping

StrepPRO™ Grouping Kit for groups A, B, C, D, F, and G

A latex agglutination test for the rapid identification of Lancefield's groups A, B, C, D, F, and G from isolated colonies of beta-hemolytic *Streptococcus* spp. and *Enterococcus faecalis*. All reagents are liquid, ready-to-use, and color coded for ease of identification.

- Greater than 98% sensitivity compared to Lancefield precipitation results. As little as one colony is needed for testing, saving time and resources without compromising accuracy.
- Complete procedure takes only two minutes

60 tests PLO30HD

Provided in kit format or separate components listed below

StrepPRO™ Extraction Reagents

Extraction reagents for the StrepPRO™ Grouping Kit. (Sold Individually)

Extraction 1, 60 tests PLO37HD

Extraction 2, 60 tests PLO38HD

Extraction 3, 60 tests PLO39HD

StrepPRO™ Grouping Reagents

(Sold Individually)

Group A, 60 tests PLO31HD

Group B, 60 tests PLO32HD

Group C, 60 tests PLO33HD

Group D, 60 tests PLO34HD

Group F, 60 tests PLO35HD

Group G, 60 tests PLO36HD

Latex Agglutination Test Cards

Disposable white cards with ten raised reaction circles used for rapid latex agglutination testing.

48/pk PLO92HD

4" Wooden Mixing Sticks

Flat-end 4" wooden mixing sticks.

75/pk PLO91HD


StrepPRO™ Polyvalent Positive Control

Positive Control for the rapid latex agglutination test for the identification of streptococcal groups.

60 tests PLO40HD

StrepPRO™ Grouping Kit for groups A, B, C, D, F, and G





HARDY DIAGNOSTICS

A Culture of Service™

40
MORE THAN  Years of
Bringing You
A Culture of Service

 **100%**
Employee Owned

**MADE IN
USA**

ESTABLISHED
1980

 Hardy Diagnostics
donates 1% of each
sale to charity.

 **bsi.** ISO
13485
Medical Devices
Quality
Management
FM 572526

Hardy Diagnostics has a Quality Management System that is certified to ISO 13485 and is a FDA licensed medical device manufacturer.

Headquarters

1430 West McCoy Lane
Santa Maria, CA 93455
800.266.2222
Sales@HardyDiagnostics.com
HardyDiagnostics.com

Distribution Centers

Santa Maria, California
Olympia, Washington
Salt Lake City, Utah
Phoenix, Arizona
Dallas, Texas
Springboro, Ohio
Lake City, Florida
Albany, New York
Raleigh, North Carolina



PREMIER®

Contract # PP-LA-640 | **Contracted Supplier**



Contract Holder

Contract # 47QSWA22D0077